


OneMarkTM

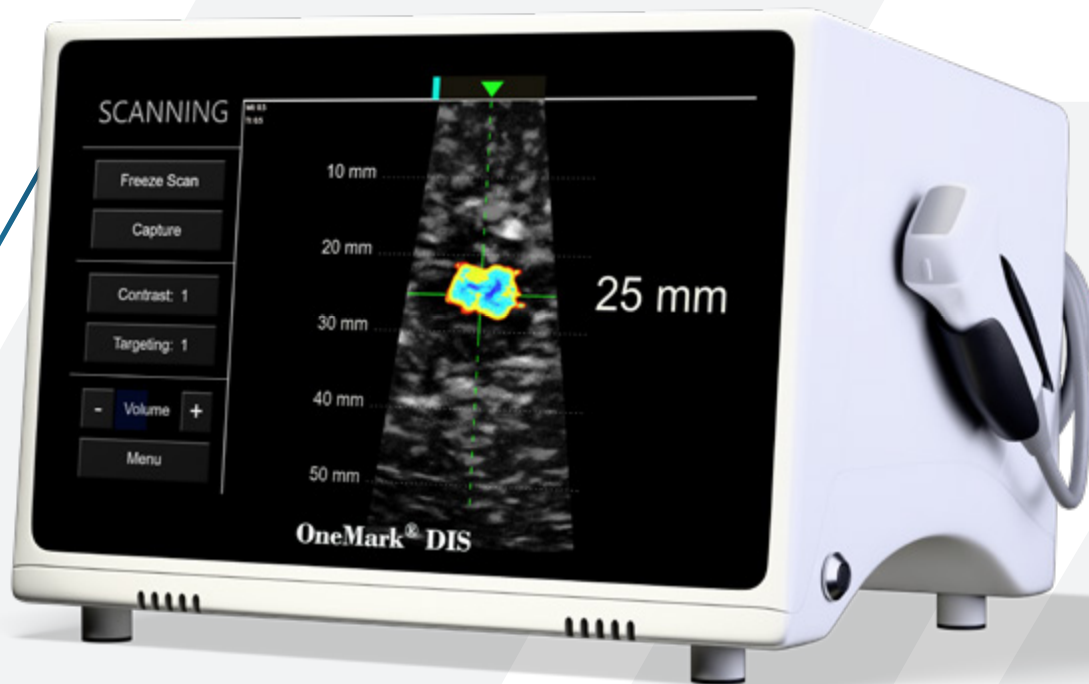

ONE PATH.
ONE LESS PROCEDURE.
ONEMARK.¹

Biopsy with OneMark

OneMark can help create a more connected path from biopsy to surgery and eliminate unnecessary procedures.¹

PLACE



Place OneMark at biopsy.

Mark biopsy sites and eliminate the need for unnecessary localization appointments.¹

NO SURGERY REQUIRED

SURGERY REQUIRED

MONITOR



Stay vigilant.

OneMark remains in place, identifying biopsy sites for further surveillance.²

SURGERY



Operate with confidence.

Progress to surgery with precision and accuracy.³

Preserve clinical options.

OneMark does not interfere with preoperative treatment, if needed, including neoadjuvant therapy and MRIs.¹

OneMark

A thoughtful design allows OneMark to remain a familiar biopsy-site marker while enabling a new kind of surgical visibility.²



Enhances Visualization

Detectable for up to 9 months, it supports biopsy through surgery, including neoadjuvant workflows.¹

Encourages Tissue Ingrowth¹

Soft and biocompatible, the marker body is made of polyethylene glycol (PEG) and collagen.

Offers Real-Time Acoustic Detection

Microspheres embedded within the marker enable system detection.¹

Supports MRI Compatibility

The OneMark marker has a clinically insignificant MRI artifact (2 mm).¹

OneMark Delivery System

To meet your clinical needs, the delivery device is available in long and short configurations ready for use at time of biopsy.¹

Position the Tip, Place the Marker

An end-deploy mechanism delivers the marker where you want it.¹

Deployment—Your Way

Our technology is designed for both single-handed use or in mammographic/stereotactic environments.²



OneMark Console


See what matters. The detection imaging console enables visualization of the marker within a tumor and surrounding tissue.¹

OneMark color overlay and crosshair indicators provide real-time visual and audio feedback to enable rapid detection and site-centering during resection.¹

Green crosshair indicators


Color overlay

Lesion and surrounding tissue

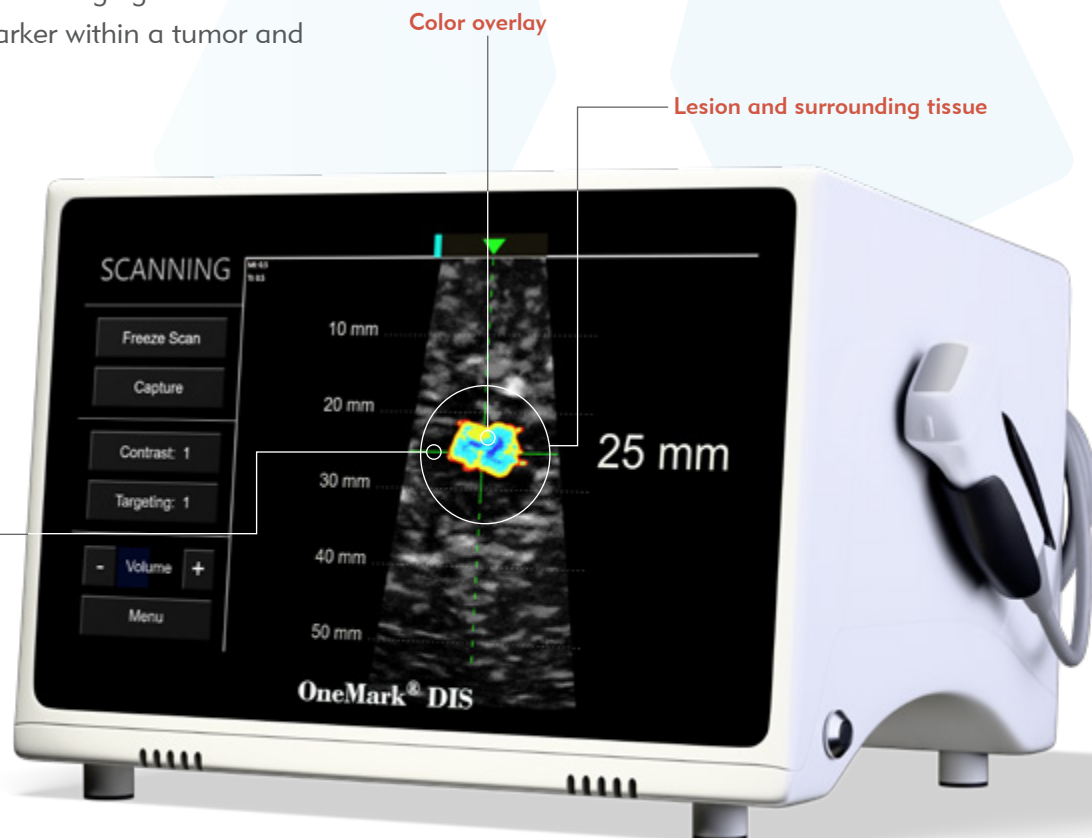

Visualize with Confidence

We are on a mission to reduce the burden cancer places on patients and their loved ones. Our comprehensive portfolio helps you localize and treat cancer with confidence, eliminating unnecessary procedures for your patients.

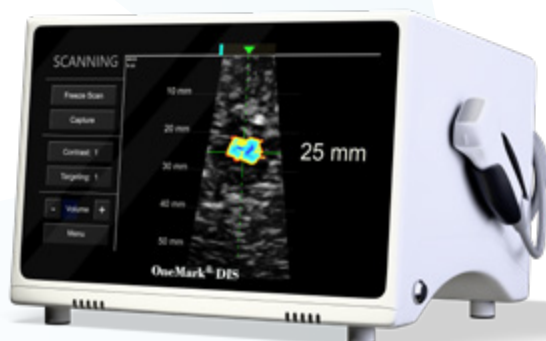


OneMark®






Contact your local Merit Medical Oncology representative to learn more.

1. Data on File.
2. OneMark Marker Instructions for Use.
3. OneMark Detection Imaging System Instructions for Use.

This brochure and other information are for the practitioner's convenience and for general information purposes only. This information does not constitute medical or legal advice, nor is it meant to endorse or guarantee the suitability of any of the referenced products or methods for any specific patient or procedure. Before using any product, refer to the Instructions for Use (IFU) for indications, contraindications, warnings, precautions, and directions for use.



Understand. Innovate. Deliver.™

meritoncology.com

Merit Medical Systems, Inc.
1600 West Merit Parkway
South Jordan, Utah 84095
1.801.253.1600
1.800.35.MERIT

Merit Medical Europe, Middle
East & Africa (EMEA)
Amerikalaan 42, 6199 AE
Maastricht-Airport
The Netherlands
+31 43 358 82 22

Merit Medical Ireland Ltd.
Parkmore Business Park West
Galway, Ireland
+353 (0) 91 703 733