**Critikit™ Model SP4500 / SP4572 In-Line Ice /
Room Temperature Injection Set
Model SP5045 Disposable In-Line Injectate Sensor Housing**

EN	INSTRUCTIONS FOR USE
FR	MODE D'EMPLOI
IT	ISTRUZIONI PER L'USO
DE	GEBRAUCHSANWEISUNG
ES	INSTRUCCIONES DE USO
PT	INSTRUÇÕES DE UTILIZAÇÃO
PTBR	INSTRUÇÕES DE USO
NL	GEbruIKSINSTRUCTIES
SV	BRUKSANVISNING
DA	BRUGSANVISNING
EL	ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ
TR	KULLANMA TALIMATLARI
PL	INSTRUKCJA UŻYTKOWANIA
CS	NÁVOD K POUŽITÍ
BG	ИНСТРУКЦИИ ЗА УПОТРЕБА
HU	HASZNÁLATI ÚTMUTATÓ
RU	РУКОВОДСТВО ПО ЭКСПЛУАТАЦИИ
LT	NAUDOJIMO INSTRUKCIJOS
LV	LIETOŠANAS PAMĀCĪBA
SK	NÁVOD NA POUŽITIE
RO	INSTRUCȚIUNI DE UTILIZARE
FI	KÄYTTÖOHJEET
HR	UPUTE ZA UPORABU
NO	BRUKSANVISNING
ET	KASUTUSJUHEND

Critikit™ Model SP4500 / SP4572 In-Line Ice / Room Temperature Injection Set

Model SP5045 Disposable In-Line Injectate Sensor Housing

INSTRUCTIONS FOR USE

DESCRIPTION

Critikit™ is a closed injectate delivery system. It connects to the CVP lumen of a thermodilution catheter to perform iced or room temperature injections to compute cardiac output.

The In-line Injectate Sensor Housing, Model SP5045 is a sensor housing designed for use with Merit or other cardiac output computers or modules, when used in conjunction with a Merit or comparable injectate sensor cable.

The Critikit™ and SP5045 when used with a Merit cardiac computer and connection cable provide measurement of injectate temperatures from 0°C to 27°C at the catheter entry site.

INTENDED USE

CRITIKIT INLINE INJECTATE DELIVERY SYSTEM

The CRITIKIT™ In-line Injectate delivery system is a closed injectate system intended for room or ice injectate temperature injection at the CVP lumen of a thermodilution catheter for computing the cardiac output.

IN-LINE SENSOR HOUSING

The in-line Sensor Housings are intended for in-line injectate temperature sensing at the catheter injection site when used with the appropriate in-line injectate sensor cables.

USER / PATIENT / CLINICAL

User: Qualified nurses, clinicians and physicians

Patient: Pediatric and adult applications

Clinical: Hospitals or appropriate clinical environments

COMPONENTS

- 1 ea. SP5045 In-Line Injectate Sensor Housing
- 1 ea. 2-way check valve assembly
- 1 ea. control syringe, 10mL
- 1 ea. IV administration spike
- 2 ea. Luer-Lok™ protective covers
- 1 ea. 3-way stopcock assembly (optional)

ADDITIONAL ITEMS FOR SETUP (SEE FIG. 1)

- SP1445 Hemopro1 Hemodynamic Profile Computer or other cardiac output computer or module for 14KΩ resistance thermodilution catheter
- SP4045 Cable or other compatible injectate sensor cable provided by Merit
- SP5555 Ice Container Bracket (Optional)
- Crushed Ice or Iced Water (Optional)

SETUP PROCEDURE

1. Use aseptic technique when handling the device.
2. Connect syringe to two-way valve assembly (Figure 3) or three-way stopcock. (Figure 2)
3. Prepare injectate solution as prescribed by a physician in a non-vented solution bag.
4. Evacuate air from injectate and IV solution bags.

NOTE: Eliminating air from IV bag will prevent air from entering the system when the solution is exhausted or when bag is inverted.
5. Close clamp on IV tubing.
6. Insert IV spike into injectate solution container.

NOTE: If a non-vented glass bottle is used, an airway needle is required to vent bottle.
7. Attach IV set to the 2-way valve assembly or three-way stopcock.
8. Verify that all components are securely connected. Hang injectate bag on IV pole. Separate small coil of tubing by gently peeling tubing apart for desired length and attachment to catheter. Place large coil in ice container.

WARNING: Unraveling the large coil may tear tubing. Tears in tubing may allow air or melted ice solution to enter the system which could be injected into the patient.
9. Open fluid control clamp on IV set.
10. Prime system by aspirating and discarding solution with syringe (see Figure 3) until entire system is free of air.

NOTE: Special care must be taken when debubbling to ensure no air is injected into patient.
11. Close fluid control clamp (and optional 3-way stopcock). Ensure that catheter CVP lumen is filled.

NOTE: A 3-way stopcock may be attached between catheter connection hub and in-line injectate sensor housing in order to add an IV solution bag for continuous flush or pressure monitoring.

WARNING: Stopcock lever must be positioned 90° for "OFF" position. Do not position stopcock at 45° to achieve an "OFF" position. 45° positioning of stopcock is imprecise and may result in contamination, bleed back or air embolism.
12. Connect injectate sensor housing and Critikit™ assembly to the prefilled CVP lumen on the PA catheter (see Figure 1). Secure injectate system to ensure minimum movement of connections.

PREPARE COOLING BATH

(Steps 13-18 are not applicable for room-temperature injectate sets)

13. Place a layer of crushed ice at bottom of container.
14. Place large coil in ice container with outlet tubing at the bottom and cover completely with ice.
15. Place unused tubing in ice container and secure tubing with slots provided at top of ice container.
16. Add ice water till full and close container lid.
17. Mount container on SP5555 Bracket and then on IV pole.
18. After proper setup and filling, injectate should be ready for use in about five minutes. Periodically check ice during use. Drain excess water into separate vessel before adding more ice.

IN-LINE SENSOR CABLE CONNECTION

19. Connect the in-line sensor cable (Model SP4045 or other compatible cable) to injectate sensor housing by pushing spring assembly into housing and tightening (see Figure 1).

20. Connect cable to Merit Hemodynamic Profile Computer, Cardiac Index Computer, Cardiac Index Module or other cardiac output computer.

NOTE: For detailed instructions refer to suppliers' instruction Manual.

CARDIAC OUTPUT MEASUREMENT

CAUTION: Catheter use must always conform to hospital policies, standards, and procedures.

21. Open roller clamp on IV set. Turn stopcock OFF to patient.

CAUTION: Use of 3mL injectate is not recommended in adults.
22. Pull back on syringe barrel to fluid-fill tubing and syringe to desired volume. Turn stopcock OFF to injectate solution.
23. Press START button on computer.
24. When you hear the tone, or when indicated by the computer, inject syringe contents rapidly into patient.
25. Read cardiac output value on computer.
26. For subsequent measurements, repeat Steps 21 through 25.
27. At conclusion of measurements, turn stopcock off to syringe and close fluid control on injectate tubing only.

WARNING

RF communicating equipment shall be maintained away from the device at a distance of 30cm (12 inches). The device shall not stack or be placed close to other communicating equipment. Electrical cables and other electronic devices may negatively affect the EMC performance.

The EMISSION characteristics of the device make it suitable for use in industrial areas and hospitals (CISPR 11 class A). Not to use in residential environment

COMPLICATIONS

Risks associated with the use of this product include: sepsis/infection, air emboli, clotted catheter and bleed-back, over infusion, and abnormal pressure readings.

For further information regarding complications, contact your Merit representative.

CAUTION: Tighten all connections before use. Do not overtighten connections as this may crack the connection leading to leaks, air embolism, bleed backs or loss of pressure waveforms.

CLINICAL BENEFITS

- Facilitate cardiac output measurement

STORAGE CONDITIONS

Store in cool dry place away from direct sunlight.

STERILE and non-pyrogenic in unopened, undamaged package. For single use only. Check integrity of the individual package before use. After use, dispose of device in a manner consistent with standard protocols for waste disposal. Do not resterilize.

Device lifetime is 72 - 96 hours base on CDC & Joint Commission Intl (JCI) recommendation.

Re-use may lead to infection or other illness/ injury.

This product contains a phthalate which may cause a risk to {children (male), pregnant or nursing women} in conjunction with prolonged duration of use.

In the EU, any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the applicable Member State.

For reordering information or assistance please contact local representative.

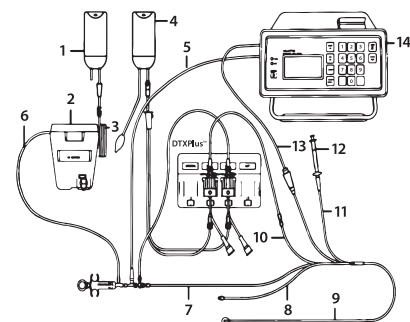
Figure 1. Closed Loop System Setup

1. Injectate solution
2. Ice Container
3. SP4500 Small Loop
4. 500 NS in a Pressure Bag
5. In-line Sensor Cable SP4045
6. Critikit™ Ice Tempkit
7. Proximal extension
8. Medication Extension
9. T.D. Criticath™ Catheter to Patient
10. Distal Extension
11. Inflation Extension
12. Volume Restricted Syringe
13. SP4035 Catheter Connector Cable
14. SP1445 Hemopro 1

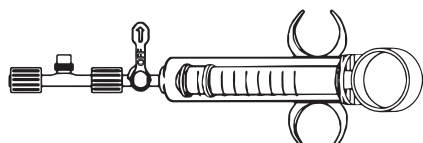
Figure 2. Sensor Housing, Stopcock and Syringe

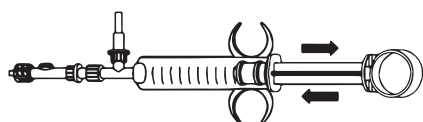
Figure 3. Filling and Debubbling with Control Syringe

	Do Not Use If Package is Damaged and Consult Instruction for Use
	Single use
	Do not resterilize
Rx ONLY	Caution: Federal (USA) law restricts this device to sale by or on the order of a physician.
	Caution
	Keep Dry
	Keep away from sunlight
	Non-pyrogenic
	Use by date: YYYY-MM-DD
	Date of Manufacture: YYYY-MM-DD
STERILE EO	Sterilized using ethylene oxide
MD	Medical Device
	Single Sterile Barrier System
UDI	Unique Device Identifier
PHT DEHP	Contains phthalates
	Warning Statement This product can expose you to chemicals including Di(2-ethylhexyl)phthalate (DEHP), which is known to the State of California to cause cancer and birth defects or other reproductive harm. For more information go to www.P65Warnings.ca.gov
	Consult Instructions for Use For electronic copy scan QR code, or go to www.merit.com/ifu and enter IFU ID. For printed copy, call U.S.A. or EU Customer Service
EC REP	Authorized Representative in European Community
	Manufacturer
REF	Catalog number
LOT	Batch code
	Defibrillation-proof type CF applied part

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Merit Medical Singapore Pte. Ltd.,
Merit Medical Ireland Ltd