DualCap® IV Pole Strips®

Disinfecting Caps for Male Luer Connectors

INSTRUCTIONS FOR USE

Each sterile DualCap® contains 70% isopropyl alcohol.

Use Aseptic Technique

DUALCAP INSTALLATION

1. Grasp the Dark Blue DualCap and remove from the IV Pole Strip®.
2. Twist the Dark Blue DualCap onto the male luer connector in a clockwise direction. There is no need to overtighten.
3. Disinfection begins in 30 seconds and continues up to 7 days under normal conditions if not removed.

TO REMOVE FOR IV LINE RECONNECTION

1. Using two hands, grasp the male luer connector and remove the Dark Blue DualCap by twisting the cap counterclockwise.
2. Reconnect the IV line.
3. Discard the DualCap.

A new Dark Blue DualCap should be placed on all male luer connectors each time the IV line is disconnected and not in use.

INDICATIONS FOR USE

When left in place for 30 seconds, the Dark Blue DualCap disinfects male luer connectors; thereafter the caps provide a physical barrier to contamination up to 7 days under normal conditions if not removed.


WARNINGS

- Do not use if package is broken, torn or punctured, or if the sponge appears to be dry.
- Do not re-use. Do not resterilize. Reuse/resterilization of DualCap may result in cross- contamination, plastic degradation, etc.
- Follow instructions for use and institutional guidelines for disposal.
- Do not store at extreme conditions.

This device is NOT made with Natural Rubber Latex.

	Caution, Consult Instructions for Use
	Caution
	Single Use. Do not reuse.
	Do not resterilize
	Use by date
	Catalog Number
	Sterilized using Irradiation
	Batch Code
	Temperature limitation
	Federal (U.S.A.) Law restricts this device to sale by or on the order of a physician
	Do not use if package is damaged


1. SEPARATE

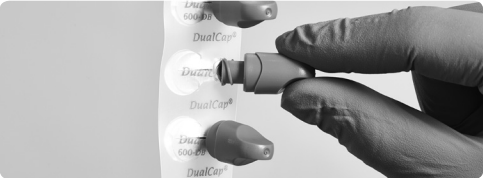2. CONNECT

3. DISINFECT

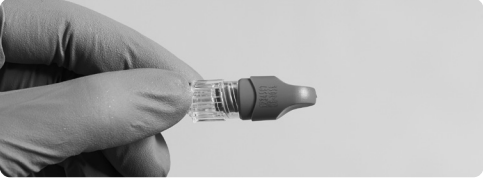Manufacturer:
Merit Medical Systems, Inc.
1600 West Merit Parkway, South Jordan, Utah
84095 U.S.A. 1-801-253-1600 U.S.A
Customer Service 1-800-356-3748

Authorized Representative:
Merit Medical Ireland Ltd,
Parkmore Business Park West, Galway, Ireland
EC Customer Service +31 43 3588222

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