

OPTIONS TO EMPOWER

Manifolds

Merit® | Devos® | TRAM®



Control, Convenience & Confidence

FIND YOUR PREFERENCE

With over 30 years' experience developing Cath Lab solutions, Merit understands that tailoring the right products to each clinician's preference is the key to optimizing outcomes and maximizing efficiency. Merit Medical offers a large selection of manifolds to meet your unique needs.

MANIFOLD FEATURES

- Clear polycarbonate body provides durability and visibility
- Ergonomic, low-torque acetyl handles designed for easy grip and smooth rotation
- Airless rotator enhances patient safety by minimizing the possibility of trapped air and creates more space for clinician's fingers
- Large bore (.093" / 2.36mm diameter) inner lumen throughout, including rotator

MANIFOLD OPTIONS

Body


- 200 & 500 PSI (14 & 34 Bar) ratings
- Right & left-handed orientation
- Block body & half body designs

Handles

- ON & OFF handle designs
- 180° & 360° rotation
- Custom colors available

Ports

- Available with up to 5 side ports
- Standard & wide port spacing
- Optional extension tubing

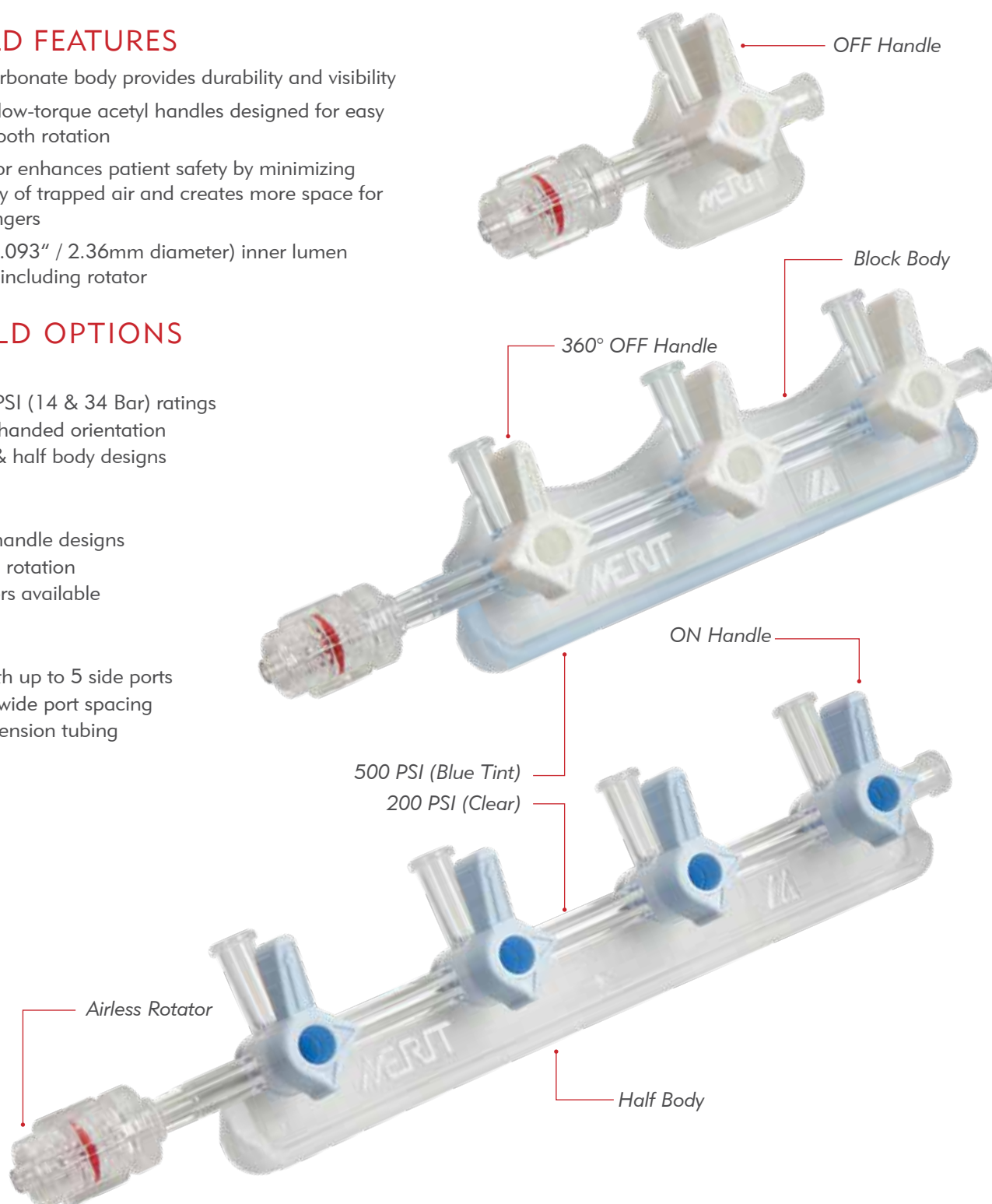



6" (15 cm) Extension Tubing



All Female Ports



Wide Port Spacing



Left-Handed Orientation



Custom Handle Colors

DEVOS™ MANIFOLD

Devos™ manifolds feature an integrated check relief valve on the port furthest from the rotator to provide directed flow control.

Integrated Check Relief Valve



TRAM® & TRAM-P™ MANIFOLDS


The TRAM® manifold with an integral transducer and 48" (122 cm) cable gives a single operator the convenience to debubble, zero, and rezero within the sterile field. Low-torque, easy-to-turn wide-spaced handles provide greater comfort, even when utilizing higher injection pressures. TRAM-P™ compensating manifolds allow pressure readings at any level.

Integrated Meritran® Pressure Transducer



ORDERING INFORMATION

			Half Body				Block Body					
			2 Port	3 Port	4 Port	1 Port	2 Port	3 Port	4 Port	5 Port		
Merit Manifolds	Right Orientation (203HF-R shown)	200 PSI (14 Bar)	OFF Handle	202HF-R	203HF-R	204HF-R	201BF-R	202BF-R	203BF-R	204BF-R		
			Wide Port	202LHWF-R	203LHWF-R	204LHWF-R						
			ON Handle	202HN-R	203HN-R	204HN-R	201BN-R	202BN-R	203BN-R	204BN-R		
			Wide Port	202LHWN-R	203LHWN-R	204LHWN-R						
		Left Orientation: Replace R with L. Not available in Wide Port. (203HF-L shown)	360° OFF Handle	202HF-RT	203HF-RT	204HF-RT	201BF-RT	202BF-RT	203BF-RT	204BF-RT		
				202HN-RT	203HN-RT	204HN-RT	201BN-RT	202BN-RT	203BN-RT	204BN-RT		
			500 PSI (34 Bar)	OFF Handle	502HF-R	503HF-R	504HF-R	501BF-R	502BF-R	503BF-R	504BF-R	
				Wide Port	502LHWF-R	503LHWF-R	504LHWF-R					
				ON Handle	502HN-R	503HN-R	504HN-R	501BN-R	502BN-R	503BN-R	504BN-R	
				Wide Port	502LHWN-R	503LHWN-R	504LHWN-R					
	360° OFF Handle	502HF-RT	503HF-RT	504HF-RT	501BF-RT	502BF-RT	503BF-RT	504BF-RT				
		502HN-RT	503HN-RT	504HN-RT	501BN-RT	502BN-RT	503BN-RT	504BN-RT				
	6" Extension Tubing (503HF-RX shown)	200 PSI (14 Bar)	OFF Handle	202HF-RX	203HF-RX	204HF-RX		202BF-RX	203BF-RX	204BF-RX		
			ON Handle	202HN-RX	203HN-RX	204HN-RX		202BN-RX	203BN-RX	204BN-RX		
			360° OFF Handle	202HF-RXT	203HF-RXT	204HF-RXT		202BF-RXT	203BF-RXT	204BF-RXT		
			360° ON Handle	202HN-RXT	203HN-RXT	204HN-RXT		202BN-RXT	203BN-RXT	204BN-RXT		
		500 PSI (34 Bar)	OFF Handle	502HF-RX	503HF-RX	504HF-RX		502BF-RX	503BF-RX	504BF-RX		
			ON Handle	502HN-RX	503HN-RX	504HN-RX		502BN-RX	503BN-RX	504BN-RX		
			360° OFF Handle	502HF-RXT	503HF-RXT	504HF-RXT		502BF-RXT	503BF-RXT	504BF-RXT		
360° ON Handle			502HN-RXT	503HN-RXT	504HN-RXT		502BN-RXT	503BN-RXT	504BN-RXT			
All Female Ports (205BN-FT shown)			200 PSI (14 Bar)	OFF Handle	202HF-F	203HF-F	204HF-F		202BF-F	203BF-F	204BF-F	205BF-F
				ON Handle	202HN-F	203HN-F	204HN-F		202BN-F	203BN-F	204BN-F	205BN-F
	360° OFF Handle	202HF-FT		203HF-FT	204HF-FT		202BF-FT	203BF-FT	204BF-FT	205BF-FT		
	360° ON Handle	202HN-FT		203HN-FT	204HN-FT		202BN-FT	203BN-FT	204BN-FT	205BN-FT		
	500 PSI (34 Bar)	OFF Handle	502HF-F	503HF-F	504HF-F		502BF-F	503BF-F	504BF-F			
		ON Handle	502HN-F	503HN-F	504HN-F		502BN-F	503BN-F	504BN-F			
		360° OFF Handle	502HF-FT	503HF-FT	504HF-FT		502BF-FT	503BF-FT	504BF-FT			
		360° ON Handle	502HN-FT	503HN-FT	504HN-FT		502BN-FT	503BN-FT	504BN-FT			
		Devos	Integral Check Relief Valve (2V3HF-RT shown)	200 PSI (14 Bar)	OFF Handle	2V2HF-R	2V3HF-R	2V4LHWF-R				
					360° OFF Handle	2V2HF-RT	2V3HF-RT	2V4LHWF-RT				
6" Extension Tubing						2V4LHWF-RXT						
TRAM	Integrated Transducer (203LTCHF-R shown)	200 PSI (14 Bar)	OFF Handle		203LTPHF-R		201LTPHF-R					
			Wide Port	202LTPHWF-R	203LTPHWF-R	204LTPHWF-R						
			Devos®			2V4LTPHWF-R						
			Devos® 360° OFF Handle			2V4LTPHWF-RT						
		500 PSI (34 Bar)	ON Handle		203LTPHN-R		201LTPHN-R					
			Wide Port	202LTPHWN-R	203LTPHWN-R	204LTPHWN-R						
			OFF Handle		503LTPHF-R		501LTPHF-R					
			Wide Port	502LTPHWF-R	503LTPHWF-R	504LTPHWF-R						
			ON Handle		503LTPHN-R		501LTPHN-R					
			Wide Port	502LTPHWN-R	503LTPHWN-R	504LTPHWN-R						
TRAM-P	Integrated Transducer with 48" Compensation Line & Clip (203LTPHN-R shown)	200 PSI (14 Bar)	OFF Handle		203LTCHF-R		201LTCHF-R					
			Wide Port	202LTCHWF-R	203LTCHWF-R	204LTCHWF-R						
			Devos®			2V4LTCHWF-R						
			Devos® 360° OFF Handle			2V4LTCHWF-RT						
		500 PSI (34 Bar)	ON Handle		203LTCHN-R		201LTCHN-R					
			Wide Port	202LTCHWN-R	203LTCHWN-R	204LTCHWN-R						
			OFF Handle		503LTCHF-R		501LTCHF-R					
			Wide Port	502LTCHWF-R	503LTCHWF-R	504LTCHWF-R						
			ON Handle		503LTCHN-R		501LTCHN-R					
			Wide Port	502LTCHWN-R	503LTCHWN-R	504LTCHWN-R						

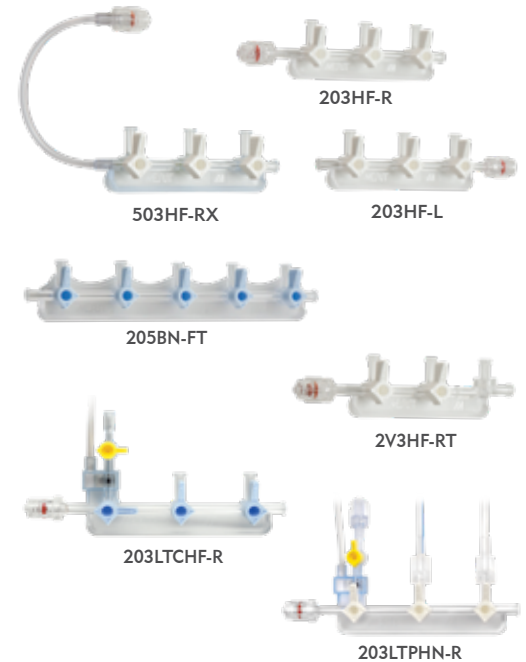

- Merit® & Devos™ Manifolds packaged 25 per box, 4 boxes per case
- TRAM® & TRAM-P™ Manifolds packaged 10 per box, 4 boxes per case

Before using, refer to Instructions for Use for indications, contraindications, warnings, precautions and directions for use.



Understand. Innovate. Deliver.™

Merit Medical Systems, Inc.
1600 West Merit Parkway
South Jordan, Utah 84095
1.801.253.1600
1.800.35.MERIT

Merit Medical Europe, Middle East, & Africa (EMEA)
Amerikalaan 42, 6199 AE Maastricht-Airport
The Netherlands
+31 43 358 82 22

Merit Medical Ireland Ltd.
Parkmore Business Park West
Galway, Ireland
+353 (0) 91 703 733

merit.com

Australia
+61 (0) 3 8373 4943

Austria
0800 295 374

Belgium
0800 72 906 (Dutch)
0800 73 172 (Français)

Brazil
+55 11 4561 - 0062

Canada
1.800.364.4370

China
+86-10 8561 0788

Denmark
80 88 00 24

Finland
0800 770 586

France
0800 91 60 30

Germany
0800 182 0871

Hong Kong
+852-2207 6888

India
+91-80-41223376

Ireland (Republic)
1800 553 163

Italy
800 897 005

Luxembourg
8002 25 22

Netherlands
0800 022 81 84

Norway
800 11629

Portugal
308 801 034

Russia
+7 495 221 89 02

Singapore
+65 67536311

South Korea
+8210 2205 0707

Spain
+34 911238406

Sweden
020 792 445

Switzerland (Deutsch)
+41 22 518 02 30

(Français)
+41 22 518 02 52

(Italiano)
+41 22 518 00 35

UK
0800 973 115