
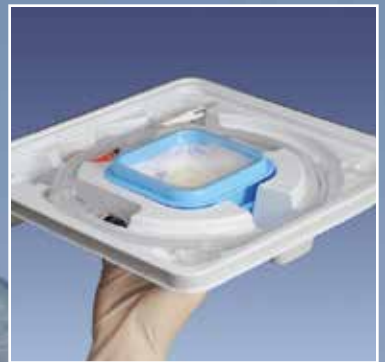


*WHEN
SECONDS
COUNT*

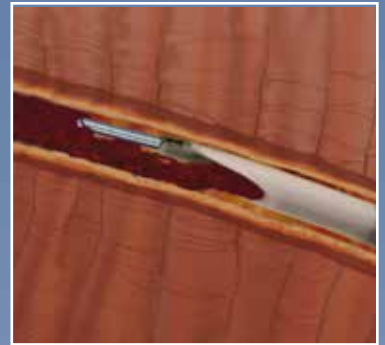


ASAP[®]
ASPIRATION CATHETER AND KIT

CUTTING EDGE CATHETER DESIGN FOR QUICK ASPIRATION

ASAP® ASPIRATION CATHETER KEY FEATURES:

- Large Aspiration Lumen to Facilitate Quick Aspiration
- Single Extrusion 100% Wire Braided Kink Resistant 140cm Catheter
- 3 Depth Positioning (nonradiopaque markers 90cm, 100cm, 110cm)
- Tapered Tip for Smooth Catheter Transition
- Hydrophilic Coating Distal 20cm of Catheter
- 12cm Rapid Exchange Length
- 6F Guide Catheter Compatibility (0.070" min. ID)
- 0.014" Guide Wire Compatibility

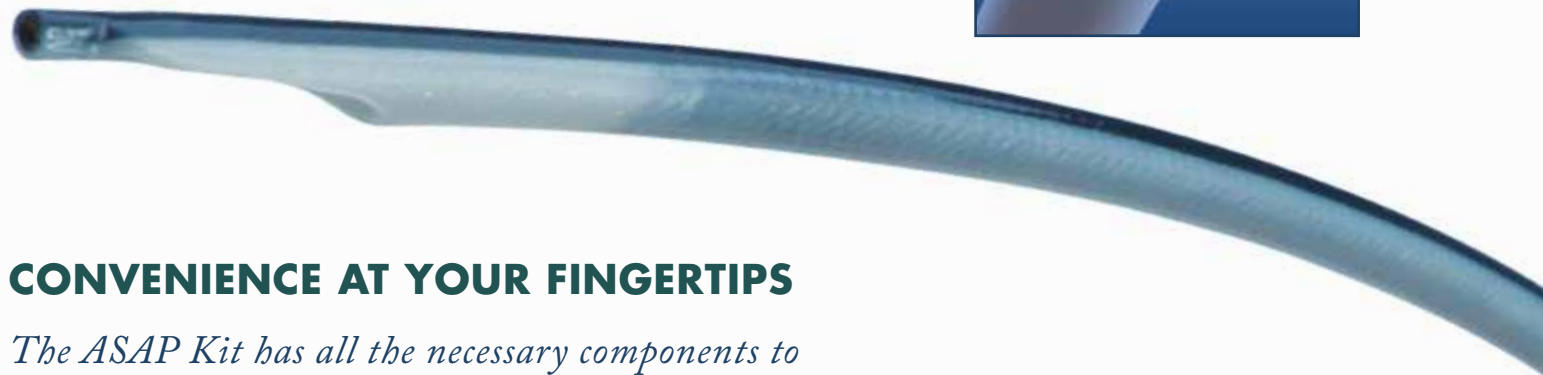

CONVENIENCE AT YOUR FINGERTIPS

The ASAP Kit has all the necessary components to perform a catheter aspiration procedure safely, quickly and effectively. With our all-in-one tray, simply remove the lid to begin the procedure.

KIT COMPONENTS:

- (1) RX 140cm Kink Resistant Aspiration Catheter
- (2) Merit VacLok® 30ml Aspiration Syringes
- (1) 4ml RXP® Flushing Syringe
- (2) Filter Baskets (70µm)
- (1) MicroStop™ Basin (Splatter proof) - 60ml fluid capacity
- (1) 21.5cm Tubing Set w/stopcock



A FAST, HANDY SOLUTION


The Merit ASAP kit is designed for the quick removal of fresh, soft emboli and thrombus from vessels of the arterial system.

THE MERIT COMPETITIVE EDGE

VENDOR	Merit ASAP	Company A	Company B
AREA OF ASPIRATION	1.56 X 10 ⁻³ sq.in.	1.08 X 10 ⁻³ sq.in.	1.44 X 10 ⁻³ sq.in.
POSITIONING MARKERS	90cm,100cm,110cm	None	95cm,105cm
FLUSH DEVICE FOR RX LUMEN	Yes - RXP	No	No
FILTER BASKETS	2 X 70 µm	2 X 40 µm	1 X 70 µm
MICROSTOP BASIN	Yes	No	No



ORDERING INFORMATION

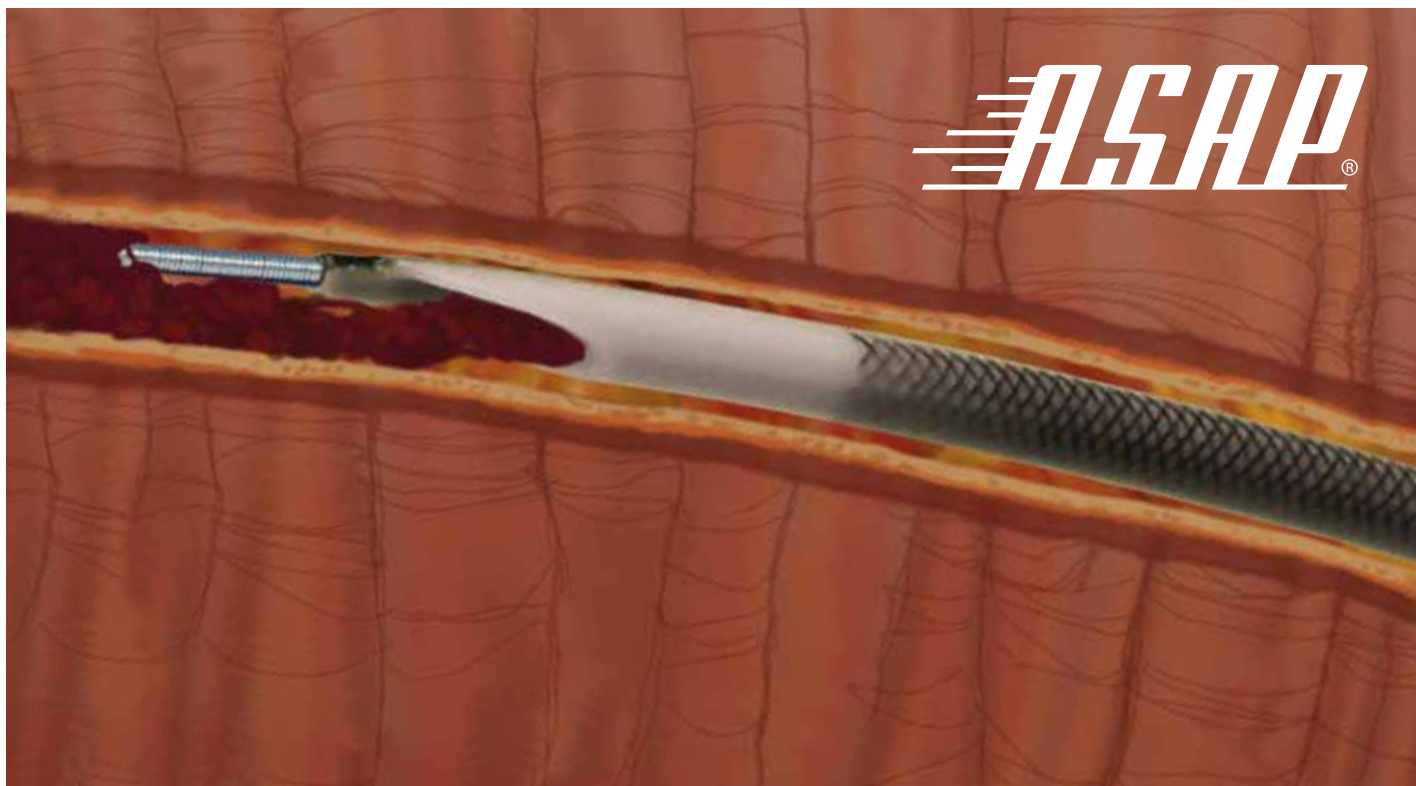

Catalog Number	Description	Box Quantity
ASAP100	ASAP Aspiration Catheter & Kit	1

U.S. Patent Nos. 6,077,698; 6,278,018; 6,603,040; 6,924,390; 7,138,541



Merit Medical Systems, Inc. • 1600 West Merit Parkway • South Jordan, Utah 84095 • 1-801-253-1600 • 1-800-35-MERIT www.merit.com

Merit Medical EUROPE, MIDDLE EAST AND AFRICA (EMEA) • Amerikalaan 42, 6199 AE Maastricht-Airport • The Netherlands • Tel: +31 43 358 82 40

Free phone for specific country:

Austria 0800 295374 • Belgium 0800 729 06 (Dutch) 0800 731 72 (French) • Denmark 80 88 00 24 • France 0800 91 60 30 • Germany 0800 182 08 71

Republic of Ireland 1800 553 163 • Luxembourg 8002 25 22 • Netherlands 0800 022 8184 • Sweden 020 79 24 45 • UK 0800 973 115 • Finland 00800 6374 8633

401951009/A ID 121312