




# Tenor<sup>®</sup>

STEERABLE GUIDEWIRE



# EmboCath<sup>®</sup> Plus

INFUSION MICROCATETER

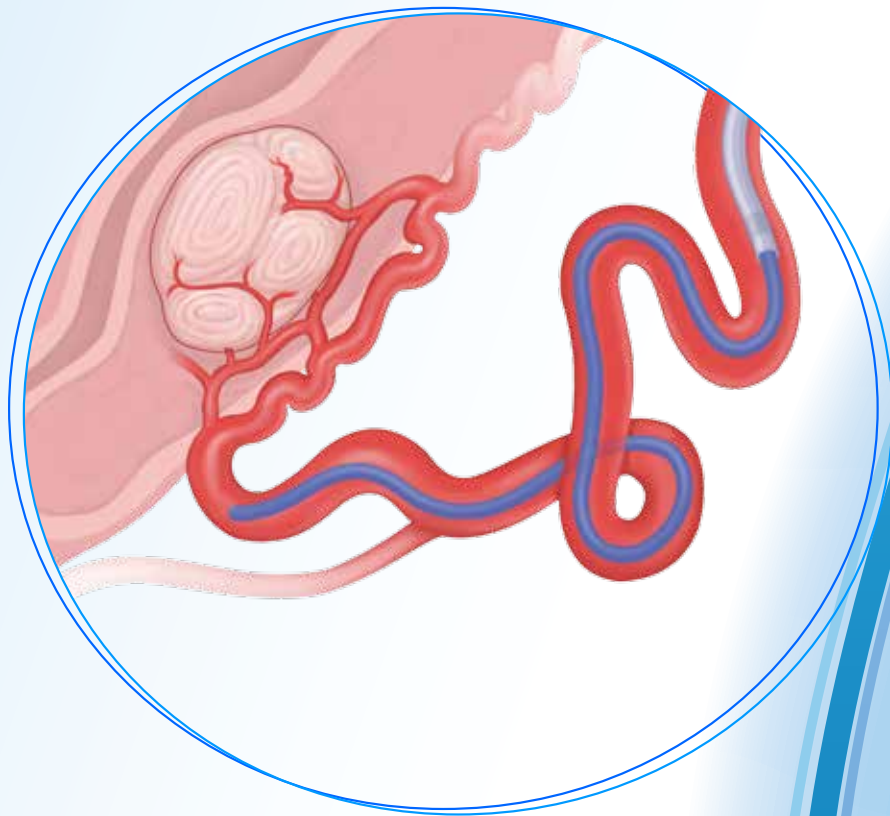

## It's All in the Delivery



# Setting the Stage for



- PFOA free PTFE on proximal segment


- **Exceptional Pushability and Torqueability** to facilitate precise navigation
- **Highly Shapeable Tip** to suit challenging anatomy
- **Radiopaque Visibility** provides for accurate placement
- **Smooth Tracking** for superb handling
- **Environmentally Friendly** - PFOA Free

- Ultra slick hydrophilic coating on distal segment


- Innovative stainless steel core wire technology

- Specially designed distal core

- Platinum coil supported by radiopaque polymer jacket

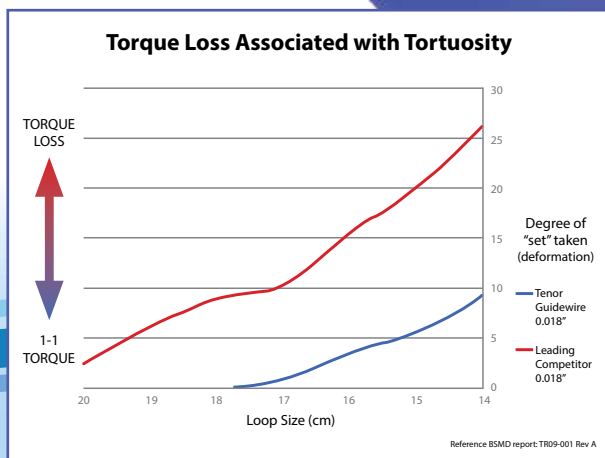

Tenor Guidewire deforms less and therefore retains torque better when manipulated into a tight loop as compared to competitors in a market evaluation.\*

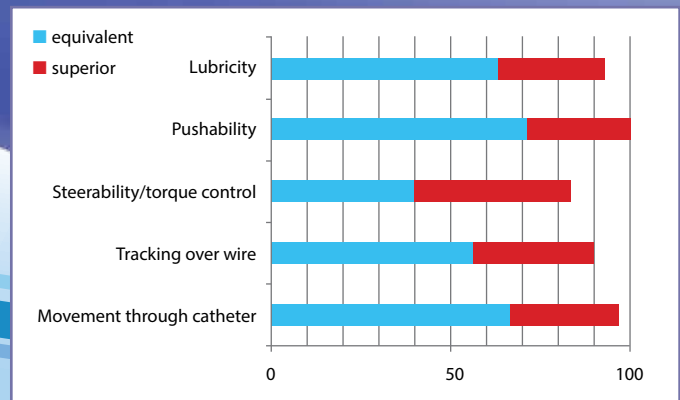

During a market evaluation more than 80% of users rated the Tenor 0.018" Guidewire as superior or equivalent to their current wire of choice across all key product attributes.\*

\*Data on file at BioSphere Medical, Inc.

# Embolotherapy

## EmboCath® Plus INFUSION MICROCATETER

### Designed Specifically for Embolization Procedures

#### • Ease of Injection

Largest inner lumen in its class  
Unique smoothly tapered hub

#### • Flexible, Kink Resistant with Excellent Pushability

Robust yet soft integrated strain relief

Pinstripe braiding for superior catheter lumen geometry

#### • Excellent Fluoroscopic Visibility

#### • Fully Chemo-compatible Material Construction




- Pinstripe® stainless steel braid technology
- Compatible with up to 700-900 µm Embosphere® Microspheres
- Hydrophilic coating
- 800psi dynamic burst pressure rating

#### • Large 0.028" Inner Diameter

- 3.0F Proximal outer diameter tapers to 2.9F
- Radiopaque marker band

### Tenor Advanced Guidewire Construction

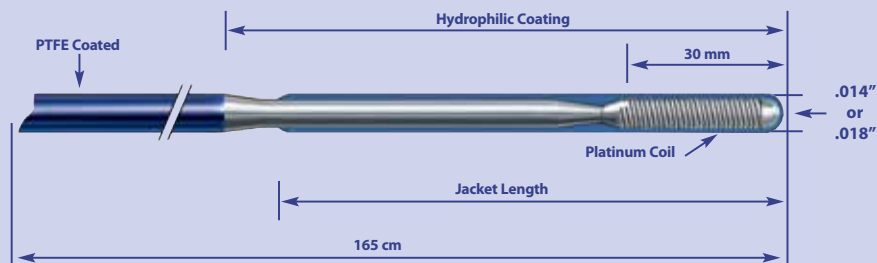

Available OD .014" and .018"  
PFOA free, PTFE on the proximal segment

## ORDERING INFORMATION

### Tenor Steerable Guidewire

Product Code	Diameter (inches)	Tip configuration	Length
TNR2411	.014"	Angled	165 cm
TNR2412	.014"	Straight	165 cm
TNR2811	.018"	Angled	165 cm
TNR2812	.018"	Straight	165 cm

### EmboCath Plus Infusion Microcatheter

Product Code	Diameter French	Length (cm)	I.D. (inches)	Tip Shape
ECP100	3.0-2.9	100	0.028	Straight
ECP135	3.0-2.9	135	0.028	Straight



Merit Medical Systems, Inc. • 1600 West Merit Parkway • South Jordan, Utah 84095 • 1-801-253-1600 • 1-800-35-MERIT www.merit.com  
 Merit Medical EUROPE, AFRICA AND THE MIDDLE EAST (EMEA) • Amerikalaan 42, 6199 AE Maastricht-Airport • The Netherlands • Tel: +31 43 358 82 22  
 BioSphere Medical, S.A. • Parc des Nations - Paris Nord 2 • 383 Rue de la Belle Etoile • 95700 Roissy en France • France  
 Free phone for specific country: Austria 0800 295 374 • Belgium 0800 72 906 (Dutch) 0800 73 172 (French) • Denmark 80 88 00 24 • France 0800 91 60 30  
 Finland 0800 770 586 • Germany 0800 182 0871 • Ireland (republic) 1800 553 163 • Italy 800 897 005 • Luxembourg 8002 25 22 • Netherlands 0800 022 81 84  
 Norway 800 11629 • Sweden 020 792 445 • UK 0800 973 115