




# PUSH TORQUE FLOW

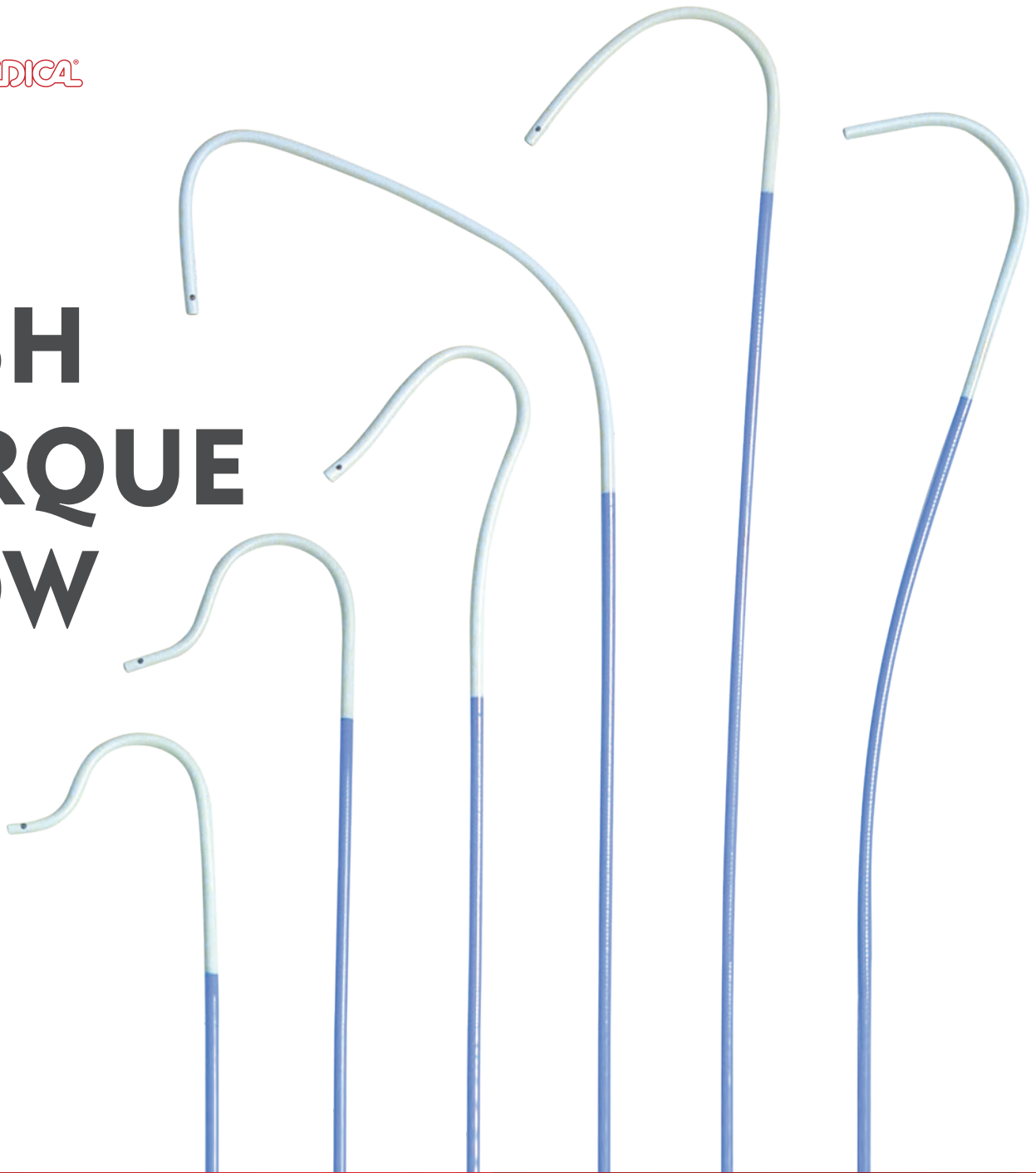

**Performa**<sup>®</sup>  
Diagnostic Peripheral Catheters

Push it to the Limit<sup>™</sup>  
with Merit's Diagnostic  
Peripheral Catheters

# Performa<sup>®</sup>

## Diagnostic Peripheral Catheters

Introducing Merit Medical's line of Diagnostic Peripheral Catheters.

Designed to provide superior torque, increased shaft strength for better pushability and large inner lumens for improved flow rates.

Now you can **PUSH IT TO THE LIMIT** in today's endovascular environment with Merit's exceptional products.

### INNER-LUMEN SIZES

- 4F .042" (1.07mm)
- 5F .046" (1.17mm)
- 6F .054" (1.37mm)
- Yields higher flow rates at lower pressures
- Superior guide wire trackability

### STRAIN RELIEF


- Increased strength at catheter/hub interface
- Color coded for easy identification

### PRESSURE

- High pressure limits at given flow rates

### WIRE BRAID DESIGN

- Greater torque for maximum control
- Improved longitudinal stability
- High degree of maneuverability
- Provides superior pushability

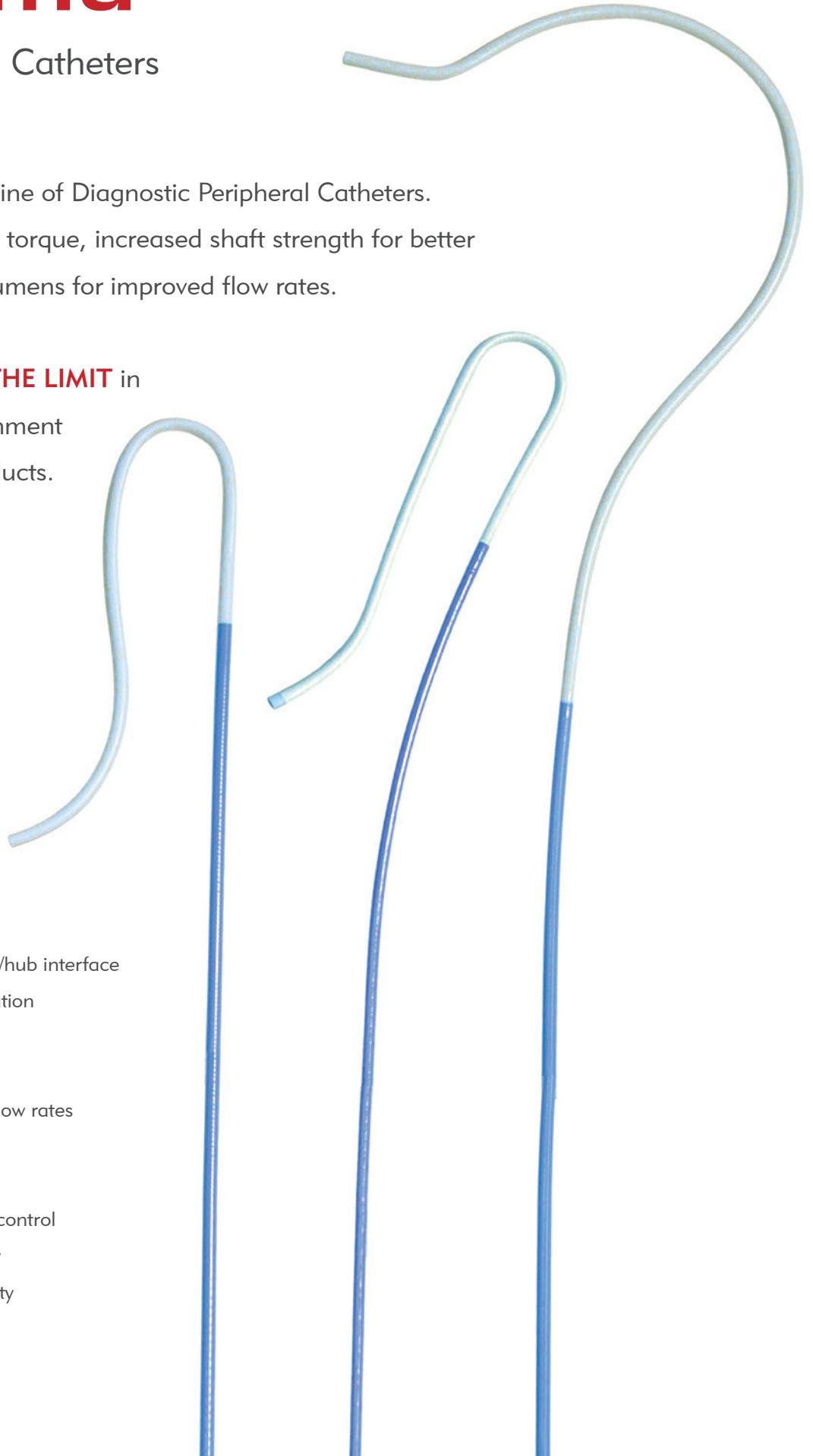

## 4F CATHETER LINE

- Accepts .038" (0.97Mm) guide wire for physician convenience
- Large inner diameter .042" (1.07Mm) provides excellent flow rates
- Bumper tip added for softness and patient safety

### WINGED POLYCARBONATE HUB

- Enhanced handling and control

### NYLON PEBAX SHAFT


- Kink Resistant
- Optimal surface smoothness
- Thermal stability at body temperature
- Thin wall design for improved flow rate

### RADIOPACITY

- Good shaft visualization under fluoroscopy
- Tip opacity ensures accurate placement and visualization

### TIP

- Soft tip for smooth atraumatic positioning
- Excellent tip memory

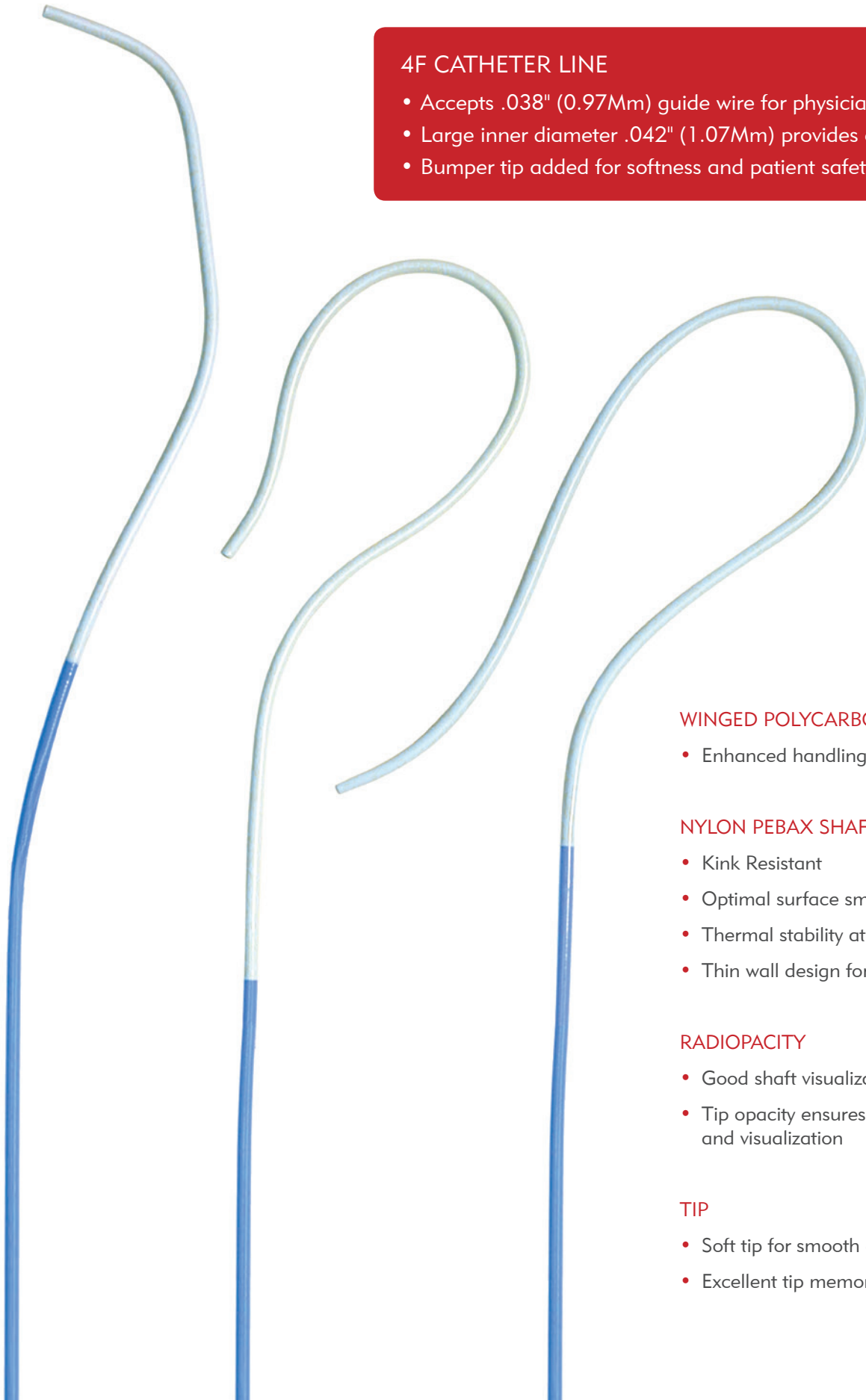

Catalog number	French Size	Internal Diameter	Length	Recommended Guide wire	# of Side Ports	Description
----------------	-------------	-------------------	--------	------------------------	-----------------	-------------

**HEADHUNTER**

7715-A0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7715-10	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5575-A3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5575-13	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	
7715-B0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7715-30	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5575-B3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5575-23	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	
7715-C0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7715-50	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5575-C3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5575-33	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	
7715-D0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7715-60	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5575-D3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5575-43	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	

**NEWTON**

7722-A0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7722-10	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5577-A3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5577-13	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	
7722-B0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7722-20	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5577-B3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5577-23	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	
7722-C0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7722-30	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5577-C3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5577-33	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	
7722-D0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7722-40	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5577-D3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5577-43	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	
7722-E0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7722-50	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5577-E3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5577-53	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	

**BENTSON**

7721-A0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7721-10	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5578-A3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5578-13	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	
7721-B0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7721-20	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5578-B3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5578-23	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	
7721-C0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7721-30	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5578-C3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5578-33	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	

1



3



5



6



1



2



3



4



5



1



2



3



Catalog number	French Size	Internal Diameter	Length	Recommended Guide wire	# of Side Ports	Description
----------------	-------------	-------------------	--------	------------------------	-----------------	-------------

**MANI**

7716-A0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7716-10	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5579-A3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5579-13	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	



MANI

**VERTEBRAL**


7717-A0	4	.042"/1.07mm	65cm	.035" (0.89mm)	-	
7717-10	4	.042"/1.07mm	65cm	.038" (0.97mm)	-	bumper tip
7717-C0	4	.042"/1.07mm	80cm	.035" (0.89mm)	-	
7717-30	4	.042"/1.07mm	80cm	.038" (0.97mm)	-	bumper tip
7717-B0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7717-20	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5583-B3	5	.046"/1.17mm	65cm	.035" (0.89mm)	-	
5583-23	5	.046"/1.17mm	65cm	.038" (0.97mm)	-	
5583-A3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5583-13	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	



VERTEBRAL

**MODIFIED CEREBRAL**

7723-A0	4	.042"/1.07mm	90cm	.035" (0.89mm)	-	
7723-10	4	.042"/1.07mm	90cm	.038" (0.97mm)	-	bumper tip
5588-A3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5588-13	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	

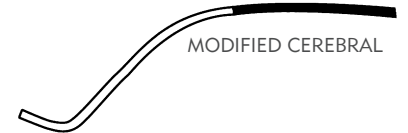

MODIFIED CEREBRAL

**BERENSTEIN**

7718-C0	4	.042"/1.07mm	40cm	.035" (0.89mm)	-	
7718-30	4	.042"/1.07mm	40cm	.038" (0.97mm)	-	bumper tip
7718-B0	4	.042"/1.07mm	65cm	.035" (0.89mm)	-	
7718-20	4	.042"/1.07mm	65cm	.038" (0.97mm)	-	bumper tip
7718-A0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7718-10	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5592-C3	5	.046"/1.17mm	40cm	.035" (0.89mm)	-	
5592-33	5	.046"/1.17mm	40cm	.038" (0.97mm)	-	
5592-B3	5	.046"/1.17mm	65cm	.035" (0.89mm)	-	
5592-23	5	.046"/1.17mm	65cm	.038" (0.97mm)	-	
5592-A3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5592-13	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	



BERENSTEIN

**STRAIGHT SELECTIVE**

7748-A0	4	.042"/1.07mm	110cm	.035" (0.89mm)	-	
7748-13	4	.042"/1.07mm	110cm	.038" (0.97mm)	-	bumper tip




STRAIGHT SELECTIVE

Catalog number	French Size	Internal Diameter	Length	Recommended Guide wire	# of Side Ports	Description
----------------	-------------	-------------------	--------	------------------------	-----------------	-------------

**MW2**

5590-13	5	.046"/1.17mm	40cm	.038" (0.97mm)	-	
5590-23	5	.046"/1.17mm	65cm	.038" (0.97mm)	-	
5590-33	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	

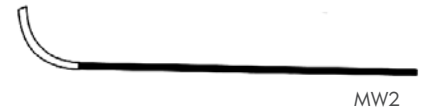

**MODIFIED MW2**

7745-A0	4	.042"/1.07mm	40cm	.035" (0.89mm)	2	
7745-10	4	.042"/1.07mm	40cm	.038" (0.97mm)	2	bumper tip
7745-B0	4	.042"/1.07mm	65cm	.035" (0.89mm)	2	
7745-C0	4	.042"/1.07mm	100cm	.035" (0.89mm)	2	



**OSBORN**

7746-A0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7746-10	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
7746-B0	4	.042"/1.07mm	65cm	.035" (0.89mm)	-	
7746-20	4	.042"/1.07mm	65cm	.038" (0.97mm)	-	bumper tip
7746-C0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7746-30	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip



**HOOK 0.8**

7711-A0	4	.042"/1.07mm	65cm	.035" (0.89mm)	2	
7711-10	4	.042"/1.07mm	65cm	.038" (0.97mm)	2	bumper tip
5585-A3	5	.046"/1.17mm	65cm	.035" (0.89mm)	2	
5585-13	5	.046"/1.17mm	65cm	.038" (0.97mm)	2	



**HOOK 1.0**

7711-B0	4	.042"/1.07mm	65cm	.035" (0.89mm)	2	
7711-20	4	.042"/1.07mm	65cm	.038" (0.97mm)	2	bumper tip
5585-B3	5	.046"/1.17mm	65cm	.035" (0.89mm)	2	
5585-23	5	.046"/1.17mm	65cm	.038" (0.97mm)	2	



**MODIFIED HOOK**

5585-C3	5	.046"/1.17mm	80cm	.035" (0.89mm)	-	
5585-D3	5	.046"/1.17mm	80cm	.035" (0.89mm)	-	
5585-D3	5	.046"/1.17mm	80cm	.035" (0.89mm)	-	

1  
2  
3



Catalog number	French Size	Internal Diameter	Length	Recommended Guide wire	# of Side Ports	Description
----------------	-------------	-------------------	--------	------------------------	-----------------	-------------

**SIMMONS**

7720-A0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7720-10	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5576-F3	5	.046"/1.17mm	65cm	.035" (0.89mm)	-	
5576-A3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5576-13	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	
7720-B0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7720-20	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5576-G3	5	.046"/1.17mm	65cm	.035" (0.89mm)	-	
5576-B3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5576-23	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	
7720-C0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7720-30	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5576-C3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5576-33	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	
7720-D0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7720-40	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5576-D3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5576-43	5	.046"/1.17mm	100cm	.038" (0.97mm)	-	

**MODIFIED SIMMONS 1.5**


7720-F0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
5576-E3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	

**MODIFIED SIMMONS**

7720-E0	4	.042"/1.07mm	80cm	.035" (0.89mm)	-	bumper tip
7720-50	4	.042"/1.07mm	80cm	.038" (0.97mm)	-	bumper tip
7720-70	4	.042"/1.07mm	110cm	.038" (0.97mm)	-	bumper tip
7720-80	4	.042"/1.07mm	120cm	.038" (0.97mm)	-	bumper tip

**COBRA 1**

7710-A0	4	.042"/1.07mm	65cm	.035" (0.89mm)	2	
7710-10	4	.042"/1.07mm	65cm	.038" (0.97mm)	2	bumper tip
7710-D0	4	.042"/1.07mm	65cm	.035" (0.89mm)	-	
7710-40	4	.042"/1.07mm	65cm	.038" (0.97mm)	-	bumper tip
5584-A3	5	.046"/1.17mm	65cm	.035" (0.89mm)	2	
5584-D3	5	.046"/1.17mm	65cm	.035" (0.89mm)	-	
5584-13	5	.046"/1.17mm	65cm	.038" (0.97mm)	-	
5584-43	5	.046"/1.17mm	65cm	.038" (0.97mm)	2	
7710-B0	4	.042"/1.07mm	65cm	.035" (0.89mm)	2	
7710-E0	4	.042"/1.07mm	65cm	.035" (0.89mm)	-	
7710-20	4	.042"/1.07mm	65cm	.038" (0.97mm)	2	bumper tip
7710-50	4	.042"/1.07mm	65cm	.038" (0.97mm)	-	bumper tip
5584-B3	5	.046"/1.17mm	65cm	.035" (0.89mm)	2	
5584-E3	5	.046"/1.17mm	65cm	.035" (0.89mm)	-	
5584-23	5	.046"/1.17mm	65cm	.038" (0.97mm)	-	
5584-53	5	.046"/1.17mm	65cm	.038" (0.97mm)	2	
7710-C0	4	.042"/1.07mm	65cm	.035" (0.89mm)	2	
7710-F0	4	.042"/1.07mm	65cm	.035" (0.89mm)	-	
7710-30	4	.042"/1.07mm	65cm	.038" (0.97mm)	2	bumper tip
7710-60	4	.042"/1.07mm	65cm	.038" (0.97mm)	-	bumper tip
5584-C3	5	.046"/1.17mm	65cm	.035" (0.89mm)	2	
5584-F3	5	.046"/1.17mm	65cm	.035" (0.89mm)	-	
5584-33	5	.046"/1.17mm	65cm	.038" (0.97mm)	-	
5584-63	5	.046"/1.17mm	65cm	.038" (0.97mm)	2	

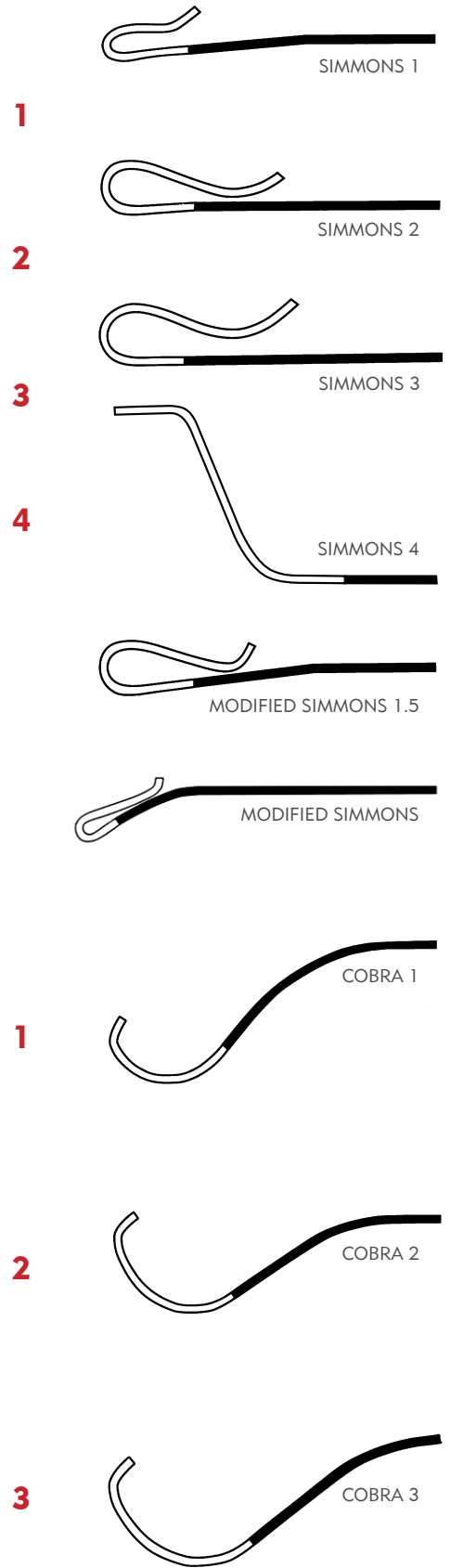

Catalog number	French Size	Internal Diameter	Length	Recommended Guide wire	# of Side Ports	Description
----------------	-------------	-------------------	--------	------------------------	-----------------	-------------

#### SHEPHERD HOOK 0.8

7712-A0	4	.042"/1.07mm	65cm	.035" (0.89mm)	2	
7712-10	4	.042"/1.07mm	65cm	.038" (0.97mm)	2	bumper tip
5586-A3	5	.046"/1.17mm	65cm	.035" (0.89mm)	2	
5586-13	5	.046"/1.17mm	65cm	.038" (0.97mm)	2	



SHEPHERD HOOK 0.8

#### SHEPHERD HOOK 1.0

7712-B0	4	.042"/1.07mm	65cm	.035" (0.89mm)	2	
7712-20	4	.042"/1.07mm	65cm	.038" (0.97mm)	2	bumper tip
5586-B3	5	.046"/1.17mm	65cm	.035" (0.89mm)	2	
5586-23	5	.046"/1.17mm	65cm	.038" (0.97mm)	2	



SHEPHERD HOOK 1.0

#### SHEPHERD HOOK 1.0 MODIFIED

7712-C0	4	.042"/1.07mm	65cm	.035" (0.89mm)	2	
7712-D0	4	.042"/1.07mm	65cm	.035" (0.89mm)	-	
7712-30	4	.042"/1.07mm	65cm	.038" (0.97mm)	2	bumper tip
7712-40	4	.042"/1.07mm	65cm	.038" (0.97mm)	-	bumper tip
5586-43	5	.046"/1.17mm	65cm	.038" (0.97mm)	2	



SHEPHERD HOOK 1.0 MODIFIED

#### RENAL DOUBLE CURVE

7713-A0	4	.042"/1.07mm	65cm	.035" (0.89mm)	2	
7713-B0	4	.042"/1.07mm	65cm	.035" (0.89mm)	-	
7713-10	4	.042"/1.07mm	65cm	.038" (0.97mm)	2	bumper tip
7713-20	4	.042"/1.07mm	65cm	.038" (0.97mm)	-	bumper tip
5587-A3	5	.046"/1.17mm	65cm	.035" (0.89mm)	2	
5587-13	5	.046"/1.17mm	65cm	.038" (0.97mm)	-	
5587-23	5	.046"/1.17mm	65cm	.038" (0.97mm)	2	



RENAL DOUBLE CURVE

#### AMIR MOTARJEME®

7747-C0	4	.042"/1.07mm	65cm	.035" (0.89mm)	-	
7747-D0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7747-E0	4	.042"/1.07mm	65cm	.035" (0.89mm)	2	
7747-F0	4	.042"/1.07mm	100cm	.035" (0.89mm)	2	
5591-A3	5	.046"/1.17mm	65cm	.035" (0.89mm)	-	
5591-B3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	
5591-C3	5	.046"/1.17mm	65cm	.035" (0.89mm)	2	
5591-D3	5	.046"/1.17mm	100cm	.035" (0.89mm)	2	



AMIR MOTARJEME

#### AMIR MOTARJEME CANE

7747-B0	4	.042"/1.07mm	65cm	.035" (0.89mm)	-	
7747-20	4	.042"/1.07mm	65cm	.038" (0.97mm)	-	bumper tip
7747-A0	4	.042"/1.07mm	100cm	.035" (0.89mm)	-	
7747-10	4	.042"/1.07mm	100cm	.038" (0.97mm)	-	bumper tip
5591-E3	5	.046"/1.17mm	65cm	.035" (0.89mm)	-	
5591-F3	5	.046"/1.17mm	100cm	.035" (0.89mm)	-	



AMIR MOTARJEME CANE



Catalog number	French Size	Internal Diameter	Length	Tip Length	Recommended Guide wire	# of Side Ports	Description
----------------	-------------	-------------------	--------	------------	------------------------	-----------------	-------------

#### HOCKEY STICK

7752-10	4	.042"/1.07mm	100cm	4.5cm	.038" (0.97mm)	-	bumper tip
7543-13	5	.046"/1.17mm	100cm	4.5cm	.038" (0.97mm)	-	
7541-13	6	.054"/1.37mm	100cm	4.5cm	.038" (0.97mm)	-	



HOCKEY STICK

Catalog number	French Size	Internal Diameter	Length	Recommended Guide wire	# of Side Ports	Description
----------------	-------------	-------------------	--------	------------------------	-----------------	-------------

#### REUTER

7714-A0	4	.042"/1.07mm	80cm	.035" (0.89mm)	2	
7714-10	4	.042"/1.07mm	80cm	.038" (0.97mm)	2	bumper tip
7714-B0	4	.042"/1.07mm	80cm	.035" (0.89mm)	-	
7714-20	4	.042"/1.07mm	80cm	.038" (0.97mm)	-	bumper tip
5589-A3	5	.046"/1.17mm	80cm	.035" (0.89mm)	-	
5589-B3	5	.046"/1.17mm	80cm	.035" (0.89mm)	2	
5589-13	5	.046"/1.17mm	80cm	.038" (0.97mm)	-	
5589-23	5	.046"/1.17mm	80cm	.038" (0.97mm)	2	



REUTER

#### MIKAELSSON

7719-A0	4	.042"/1.07mm	80cm	.035" (0.89mm)	-	
7719-10	4	.042"/1.07mm	80cm	.038" (0.97mm)	-	bumper tip
5594-13	5	.046"/1.17mm	80cm	.038" (0.97mm)	-	



MIKAELSSON

#### KA CATHETERS

5582-23	5	.046"/1.17mm	65cm	.038" (0.97mm)	-	
5582-33	5	.052"/1.32mm	65cm	.038" (0.97mm)	-	high flow non-braided



KA CATHETERS

#### KA2 CATHETERS

5582-13	5	.046"/1.17mm	40cm	.038" (0.97mm)	-	
5582-43	5	.046"/1.17mm	65cm	.038" (0.97mm)	-	
7514-13	6	.054"/1.37mm	40cm	.038" (0.97mm)	-	
7514-23	6	.054"/1.37mm	65cm	.038" (0.97mm)	-	



KA2 CATHETERS

#### DVS A1 (10cm taper tip)

5595-13	5	.046"/1.12mm	65cm	.035" (0.89mm)	-	
5595-33	5	.046"/1.12mm	100cm	.035" (0.89mm)	-	
5595-23	5	.046"/1.12mm	65cm	.038" (0.97mm)	-	
5595-43	5	.046"/1.12mm	100cm	.038" (0.97mm)	-	



DVS A1

#### DVS A2 (15cm taper tip)

5595-A3	5	.046"/1.12mm	65cm	.035" (0.89mm)	-	
5595-C3	5	.046"/1.12mm	100cm	.035" (0.89mm)	-	
5595-B3	5	.046"/1.12mm	65cm	.038" (0.97mm)	-	
5595-D3	5	.046"/1.12mm	100cm	.038" (0.97mm)	-	



DVS A2

Catalog number	French Size	Internal Diameter	Length	Recommended Guide wire	# of Side Ports	Description
----------------	-------------	-------------------	--------	------------------------	-----------------	-------------

#### UHF SHEPHERD FLUSH

7602-A1	5	.052"/1.32mm	65cm	.035" (0.89mm)	6	non-braided
7602-B1	5	.052"/1.32mm	90cm	.035" (0.89mm)	6	non-braided
7602-C1	5	.052"/1.32mm	65cm	.038" (0.97mm)	6	non-braided
7602-D1	5	.052"/1.32mm	90cm	.038" (0.97mm)	6	non-braided



#### ULTRA BOLUS FLUSH

7674-11	5	.052"/1.32mm	65cm	.035" (0.89mm)	10	non-braided
7674-21	5	.052"/1.32mm	90cm	.035" (0.89mm)	10	non-braided
7674-31	5	.052"/1.32mm	65cm	.038" (0.97mm)	10	non-braided
7674-41	5	.052"/1.32mm	90cm	.038" (0.97mm)	10	non-braided



#### PIGTAIL FLUSH

7749-A0	4	.042"/1.07mm	80cm	.035" (0.89mm)	10	braided
7749-40	4	.042"/1.07mm	80cm	.038" (0.97mm)	10	braided



#### ULTRA HIGH FLOW PIGTAIL

7602-41	5	.052"/1.32mm	65cm	.035" (0.89mm)	10	non-braided
7602-51	5	.052"/1.32mm	90cm	.035" (0.89mm)	10	non-braided
7602-61	5	.052"/1.32mm	110cm	.035" (0.89mm)	10	non-braided
7602-71	5	.052"/1.32mm	65cm	.038" (0.97mm)	10	non-braided
7602-81	5	.052"/1.32mm	90cm	.038" (0.97mm)	10	non-braided
7602-91	5	.052"/1.32mm	110cm	.038" (0.97mm)	10	non-braided



#### STRAIGHT FLUSH

7673-11	5	.052"/1.32mm	65cm	.035" (0.89mm)	10	non-braided
7673-21	5	.052"/1.32mm	90cm	.035" (0.89mm)	10	non-braided
7673-31	5	.052"/1.32mm	65cm	.035" (0.89mm)	-	non-braided
7673-41	5	.052"/1.32mm	65cm	.038" (0.97mm)	10	non-braided
7673-51	5	.052"/1.32mm	90cm	.038" (0.97mm)	10	non-braided
7673-61	5	.052"/1.32mm	65cm	.038" (0.97mm)	-	non-braided



#### MODIFIED HOOK FLUSH

7603-A1	5	.052"/1.32mm	65cm	.035" (0.89mm)	6	non-braided
7603-B1	5	.052"/1.32mm	90cm	.035" (0.89mm)	6	non-braided
7603-C1	5	.052"/1.32mm	65cm	.038" (0.97mm)	6	non-braided
7603-D1	5	.052"/1.32mm	90cm	.038" (0.97mm)	6	non-braided



Catalog number	French Size	Internal Diameter	Length	Recommended Guide wire	# of Side Ports	Bands/Spacing	Tip Style
----------------	-------------	-------------------	--------	------------------------	-----------------	---------------	-----------

**PERFORMA VESSEL SIZING All side ports between tip and first band**

7602-2M	5	.052"/1.32mm	100cm	.035" (0.89mm)	10	2 bands 2cm spacing	Pigtail
7602-2M65	5	.052"/1.32mm	65cm	.035" (0.89mm)	10	2 bands 2cm spacing	Pigtail
7674-2M	5	.052"/1.32mm	100cm	.035" (0.89mm)	10	2 bands 2cm spacing	Ultra Bolus
7674-2M65	5	.052"/1.32mm	65cm	.035" (0.89mm)	10	2 bands 2cm spacing	Ultra Bolus
7673-2M	5	.052"/1.32mm	100cm	.035" (0.89mm)	10	2 bands 2cm spacing	Straight
7673-2M65	5	.052"/1.32mm	65cm	.035" (0.89mm)	10	2 bands 2cm spacing	Straight
7603-2M100	5	.052"/1.32mm	100cm	.035" (0.89mm)	6	2 bands 2cm spacing	Mod. Hook
7603-2M65	5	.052"/1.32mm	65cm	.035" (0.89mm)	6	2 bands 2cm spacing	Mod. Hook
7742-2M65SH	4	.041"/1.04mm	65cm	.035" (0.89mm)	5	2 (2cm apart)	Pigtail
7742-2M100SH	4	.041"/1.04mm	100cm	.035" (0.89mm)	5	2 (2cm apart)	Pigtail
7744-2M65SH	4	.041"/1.04mm	65cm	.035" (0.89mm)	5	2 (2cm apart)	Straight
7744-2M100SH	4	.041"/1.04mm	100cm	.035" (0.89mm)	5	2 (2cm apart)	Straight
7743-2M65SH	4	.041"/1.04mm	65cm	.035" (0.89mm)	5	2 (2cm apart)	Ultra Bolus
7743-2M100SH	4	.041"/1.04mm	100cm	.035" (0.89mm)	5	2 (2cm apart)	Ultra Bolus
7755-2M65	4	.041"/1.04mm	65cm	.035" (0.89mm)	6	2 (2cm apart)	Mod. Hook
7755-2M100	4	.041"/1.04mm	100cm	.035" (0.89mm)	6	2 (2cm apart)	Mod. Hook
7602-11M	5	.052"/1.32mm	100cm	.035" (0.89mm)	10	11 bands 1cm spacing	Pigtail
7602-11M65	5	.052"/1.32mm	65cm	.035" (0.89mm)	10	11 bands 1cm spacing	Pigtail
7673-11M100SH	5	.052"/1.32mm	100cm	.035" (0.89mm)	10	11 bands 1cm spacing*	Straight
7673-11M65SH	5	.052"/1.32mm	65cm	.035" (0.89mm)	10	11 bands 1cm spacing*	Straight
7603-11M65	5	.052"/1.32mm	65cm	.035" (0.89mm)	6	11 bands 1cm spacing	Mod. Hook
7603-11M100	5	.052"/1.32mm	100cm	.035" (0.89mm)	6	11 bands 1cm spacing	Mod. Hook
7742-20M65SH	4	.041"/1.04mm	65cm	.035" (0.89mm)	5	20 bands 1cm spacing	Pigtail
7742-20M100SH	4	.041"/1.04mm	100cm	.035" (0.89mm)	5	20 bands 1cm spacing	Pigtail
7602-20M	5	.052"/1.32mm	100cm	.035" (0.89mm)	10	20 bands 1cm spacing	Pigtail
7602-20M65	5	.052"/1.32mm	65cm	.035" (0.89mm)	10	20 bands 1cm spacing	Pigtail
7602-20M65SH	5	.052"/1.32mm	65cm	.035" (0.89mm)	5	20 bands 1cm spacing*	Pigtail
7602-20M100SH	5	.052"/1.32mm	100cm	.035" (0.89mm)	5	20 bands 1cm spacing*	Pigtail
7673-20M65SH	5	.052"/1.32mm	65cm	.035" (0.89mm)	10	20 bands 1cm spacing*	Straight
7673-20M100SH	5	.052"/1.32mm	100cm	.035" (0.89mm)	10	20 bands 1cm spacing*	Straight
7603-20M100SH	5	.052"/1.32mm	100cm	.035" (0.89mm)	6	20 bands 1cm spacing*	Mod. Hook
7603-20M65SH	5	.052"/1.32mm	65cm	.035" (0.89mm)	6	20 bands 1cm spacing*	Mod. Hook



STRAIGHT

**PERFORMA VESSEL SIZING All side ports between bands**

56535AVS2PIG	5	.052"/1.32mm	65cm	.035" (0.89mm)	10	2 (2cm apart)	Pigtail
510035AVS2PIG	5	.052"/1.32mm	100cm	.035" (0.89mm)	10	2 (2cm apart)	Pigtail
56535AVS2STR	5	.052"/1.32mm	65cm	.035" (0.89mm)	10	2 (2cm apart)	Straight
510035AVS2STR	5	.052"/1.32mm	100cm	.035" (0.89mm)	10	2 (2cm apart)	Straight
56535AVS2MHK	5	.052"/1.32mm	65cm	.035" (0.89mm)	6	2 (2cm apart)	Mod Hook
510035AVS2MHK	5	.052"/1.32mm	100cm	.035" (0.89mm)	6	2 (2cm apart)	Mod Hook
46535AVS20MHK	4	.041"/1.04mm	65cm	.035" (0.89mm)	6	20 (1cm apart)	Mod Hook
56535AVS20PIG	5	.052"/1.32mm	65cm	.035" (0.89mm)	10	20 (1cm apart)	Pigtail
510035AVS20PIG	5	.052"/1.32mm	100cm	.035" (0.89mm)	10	20 (1cm apart)	Pigtail
56535AVS20MHK	5	.052"/1.32mm	65cm	.035" (0.89mm)	6	20 (1cm apart)	Mod Hook
510035AVS20MHK	5	.052"/1.32mm	100cm	.035" (0.89mm)	6	20 (1cm apart)	Mod Hook



**PERFORMA PEDIATRIC VESSEL SIZING**

7757-2M65	4	.041"/1.04mm	65cm	.035" (0.89mm)	5	2 bands (1cm apart)	Pigtail
7757-2M80	4	.041"/1.04mm	80cm	.035" (0.89mm)	5	2 bands (1cm apart)	Pigtail
7757-2M65SH	4	.041"/1.04mm	65cm	.035" (0.89mm)	5	2 bands (1cm apart)*	Pigtail
7757-2M80SH	4	.041"/1.04mm	80cm	.035" (0.89mm)	5	2 bands (1cm apart)*	Pigtail

\* Bands are 16mm closer to pigtail curve

UNDERSTAND.  
INNOVATE.  
DELIVER.™

Before using refer to Instructions for Use for indications, contraindications, warnings, precautions, and directions for use.



Understand. Innovate. Deliver.™

Merit Medical Systems, Inc.  
1600 West Merit Parkway  
South Jordan, Utah 84095  
1.801.253.1600  
1.800.35.MERIT

Merit Medical Europe, Middle East, & Africa (EMEA)  
Amerikalaan 42, 6199 AE Maastricht-Airport  
The Netherlands  
+31 43 358 82 22

Merit Medical Ireland Ltd.  
Parkmore Business Park West  
Galway, Ireland  
+353 (0) 91 703 733

[merit.com](http://merit.com)

Australia  
+61 (0) 3 8373 4943

Austria  
0800 295 374

Belgium  
0800 72 906 (Dutch)  
0800 73 172 (Français)

Brazil  
+55 11 4561 - 0062

Canada  
1.800.364.4370

China  
+86-10 8561 0788

Denmark  
80 88 00 24

Finland  
0800 770 586

France  
0800 91 60 30

Germany  
0800 182 0871

Hong Kong  
+852-2207 6888

India  
+91-80-41223376

Ireland (Republic)  
1800 553 163

Italy  
800 897 005

Luxembourg  
8002 25 22

Netherlands  
0800 022 81 84

Norway  
800 11629

Portugal  
308 801 034

Russia  
+7 495 221 89 02

Singapore  
+65 67536311

South Korea  
+8210 2205 0707

Spain  
+34 911238406

Sweden  
020 792 445

Switzerland  
(Deutsch)  
+41 22 518 02 30

(Français)  
+41 22 518 02 52

(Italiano)  
+41 22 518 00 35

UK  
0800 973 115