



PERFORMANCE AT EVERY CURVE

DYNAMIS ΔV ™

High Pressure PTA Dilatation Catheters

Balloons that conform to the
central venous anatomy

DYNAMIS Δ V™

BUILT FOR ARTERIOVENOUS DIALYSIS FISTULAE

The majority of dialysis access-related blockages are located in the venous system, where the anatomy is tortuous and complex. The Dynamis AV™ facilitates the opening of blockages located in the arteriovenous system of dialysis patients. It is built on an **0.035" OVER-THE-WIRE PLATFORM** with up to **27 atm RBP** and a **75 cm CATHETER SHAFT LENGTH**.*

Rated Burst Pressure

14 – 27 atm to dilate resistant lesions

Radiopaque Marker Bands

Facilitate precise balloon positioning


Semi-compliant Nylon Balloon

Promotes vessel conformability and resists vessel straightening

Vessel wall


*Indicated for use in percutaneous transluminal angioplasty of the femoral, iliac and renal arteries, and for the treatment of obstructive lesions of native or synthetic arteriovenous dialysis fistulae.

Actual photograph enlarged to show detail

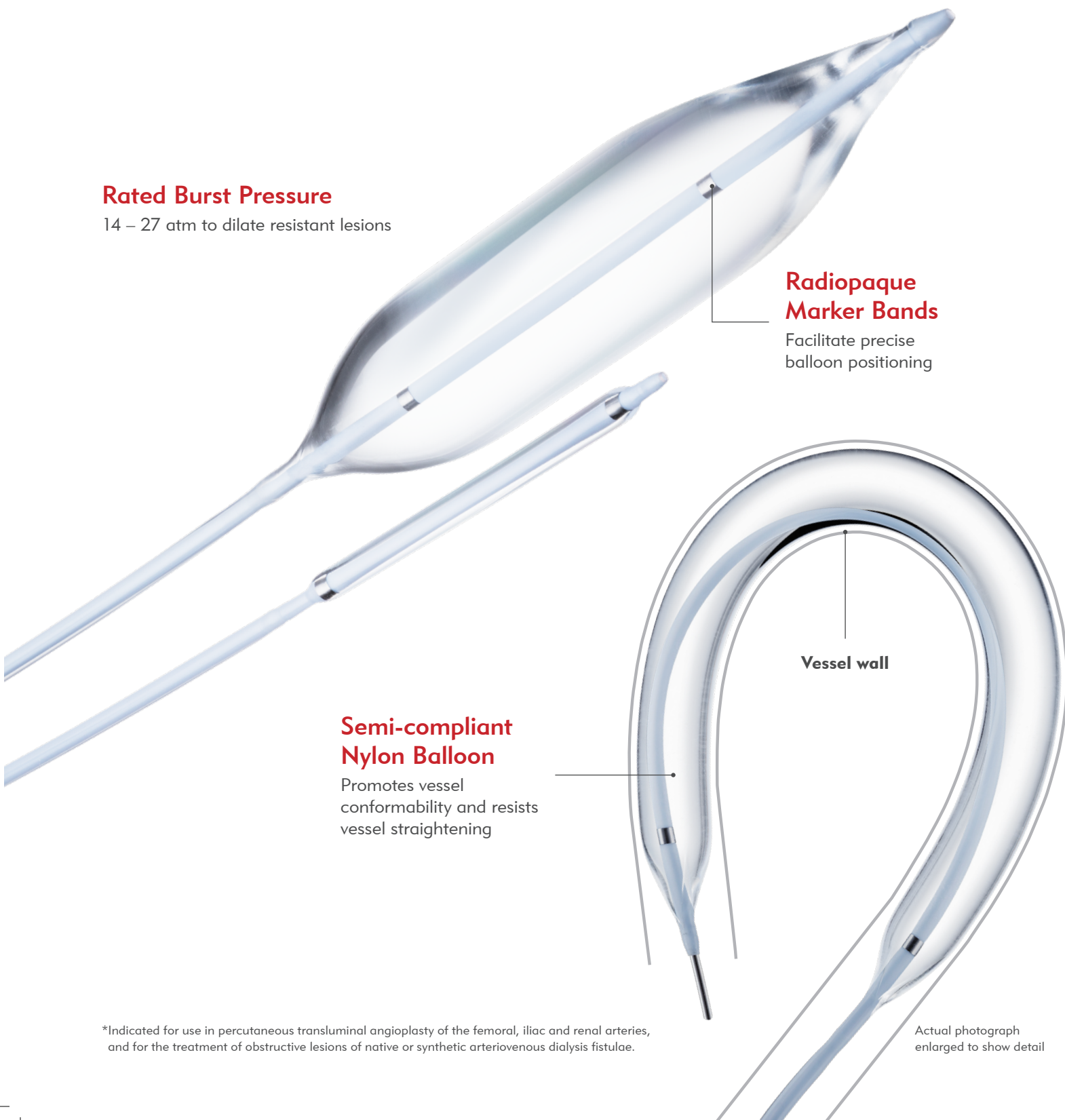


Broad Balloon Diameter Range

3 mm – 16 mm

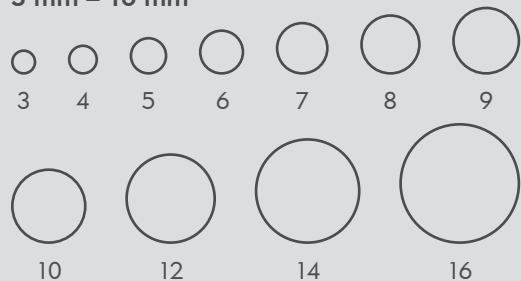

Coaxial Catheter Shaft

Supports rapid inflation, deflation, trackability, and pushability

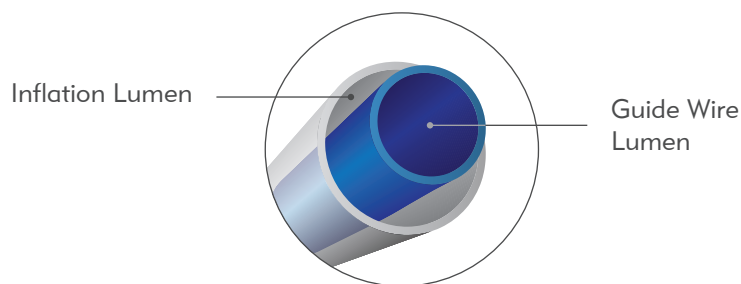

Pre-installed Re-wrap Tool

Eases reinsertion through the sheath




The high pressure Dynamis AV™ balloon is designed to dilate resistant lesions while conforming to the vessel, as demonstrated in the case studies below:**


Severe Juxta-anastomotic Stenosis

- 66 Male with ESRD, CAD, COPD
- L RC AVF with difficult cannulation, thrombus present
- Exam—poor augmentation


Dynamis AV 7x40 22 atm

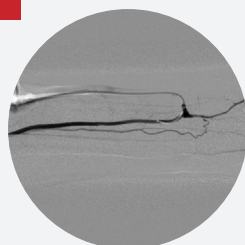

Before

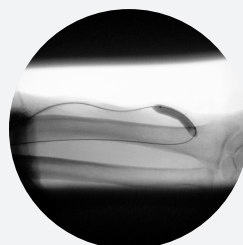

During

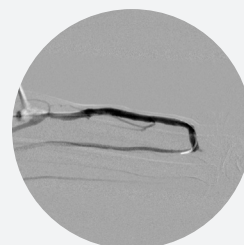


After


99% Cephalic Arch Stenosis

- 62 Female with ESRD secondary to FSGS
- L BC AVF
- Presents with prolonged bleeding
- Exam reveals a pulsatile fistula that does not collapse with arm elevation


Dynamis AV 12x40 18 atm

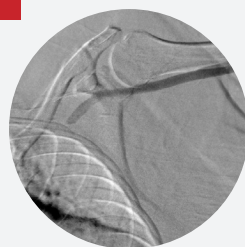

Before

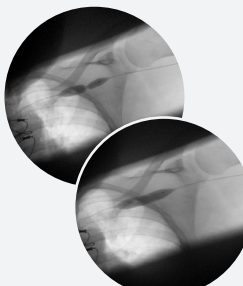

During

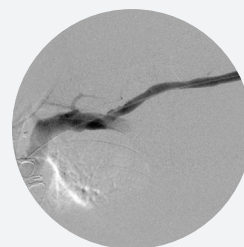

After

**Used by permission



Dynamis AV™ is part of the Merit family of Dialysis Access Products: introducers, guide wires, inflation devices, snare systems, suture and compression devices.

Please visit [Merit.com/Dialysis](https://www.merit.com/Dialysis) or contact your Merit Sales Professional to learn more.



DYNAMIS ΔV ™

High Pressure PTA Dilatation Catheters

75 CM CATHETER SHAFT LENGTH, .035" GUIDEWIRE COMPATIBILITY

Model Number	Balloon Diameter (mm)	Balloon Length (mm)	Nominal Pressure (atm)	Rated Burst (atm)	Sheath Size (F)
D0304062775	3	40	14	27	6
D0402062775	4	20	14	27	6
D0404062775	4	40	14	27	6
D0502062775	5	20	14	27	6
D0504062775	5	40	14	27	6
D0602062575	6	20	14	25	6
D0604062575	6	40	14	25	6
D0610062575	6	100	14	25	6
D0702062575	7	20	14	25	6
D0704062575	7	40	14	25	6
D0704062375S	7	40	14	23	6
D0706062375	7	60	14	23	6
D0710062375	7	100	14	23	6
D0802062475	8	20	14	24	6
D0804062475	8	40	14	24	6
D0804062275S	8	40	14	22	6
D0806062275	8	60	14	22	6
D0808062275	8	80	14	22	6
D0904072075	9	40	12	20	7
D1004072075	10	40	12	20	7
D1204071875	12	40	12	18	7
D1404081675	14	40	10	16	8
D1604081475	16	40	10	14	8

Balloon Length

20 mm – 100 mm

Before using refer to Instructions for Use for indications, contraindications, warnings, precautions, and directions for use.



Understand. Innovate. Deliver.™

Merit Medical Systems, Inc.
1600 West Merit Parkway
South Jordan, Utah 84095
1.801.253.1600
1.800.35.MERIT

merit.com

Merit Medical Europe, Middle East, & Africa (EMEA)
Amerikalaan 42, 6199 AE Maastricht-Airport
The Netherlands
+31 43 358 82 22

Merit Medical Ireland Ltd.
Parkmore Business Park West
Galway, Ireland
+353 (0) 91 703 733

Australia
+61 (0) 3 8373 4943

Austria
0800 295 374

Belgium
0800 72 906 (Dutch)
0800 73 172 (Français)

Brazil
+55 11 4561 - 0062

Canada
1.800.364.4370

China
+86-10 8561 0788

Denmark
80 88 00 24

Finland
0800 770 586

France
0800 91 60 30

Germany
0800 182 0871

Hong Kong
+852-2207 6888

India
+91-80-41223376

Ireland (Republic)
1800 553 163

Italy
800 897 005

Luxembourg
8002 25 22

Netherlands
0800 022 81 84

Norway
800 11629

Portugal
308 801 034

Russia
+7 495 221 89 02

Singapore
+65 67536311

South Korea
+8210 2205 0707

Spain
+34 911238406

Sweden
020 792 445

Switzerland
(Deutsch)
+41 22 518 02 30

(Français)
+41 22 518 02 52

(Italiano)
+41 22 518 00 35

UK
0800 973 115