



ADVANCING OPTIONS FOR RIGHT SIDED ACCESS

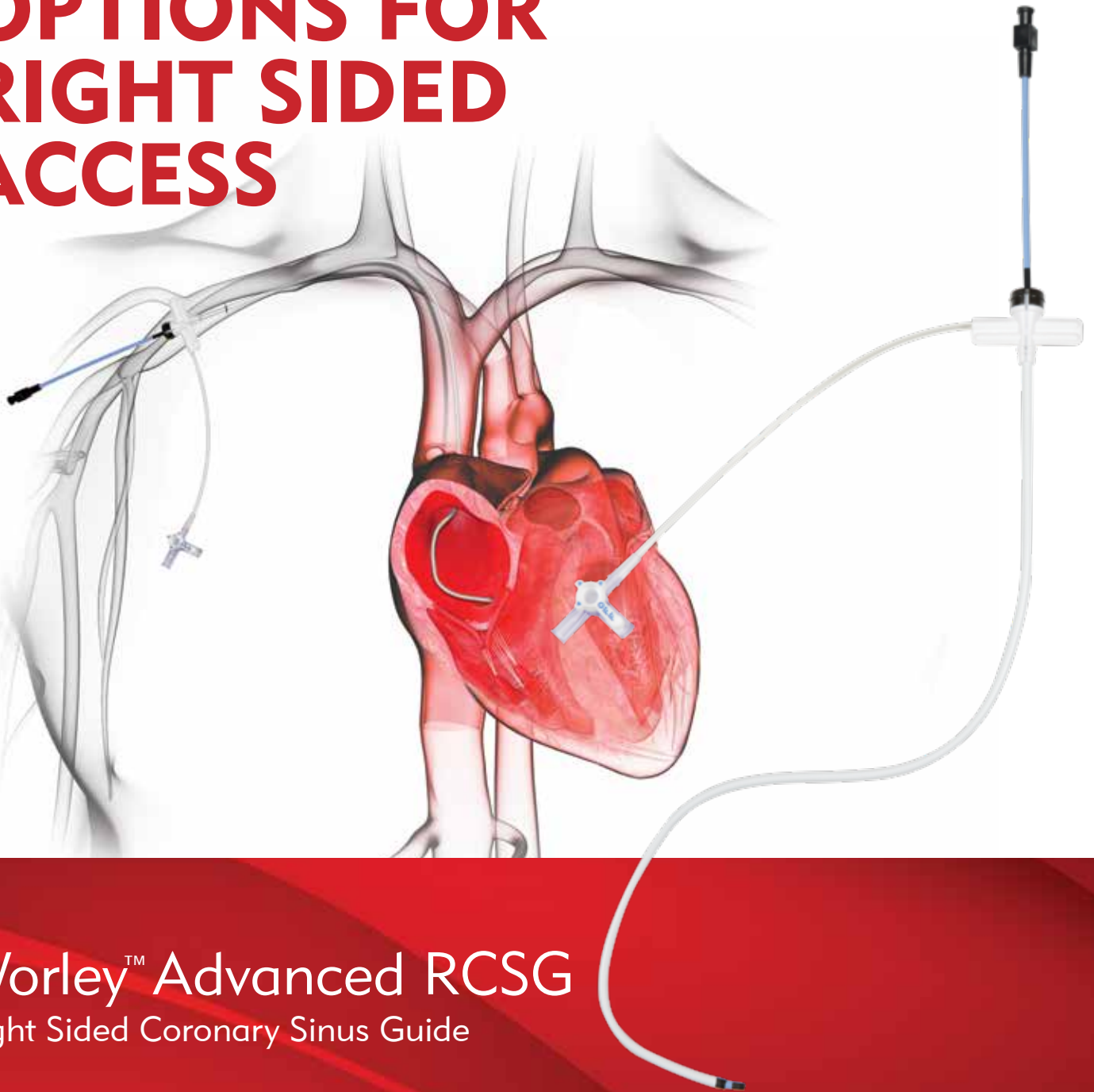

Worley™ Advanced RCSG

Right Sided Coronary Sinus Guide

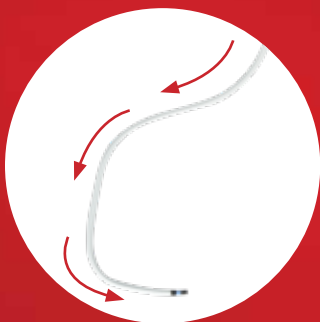

Unique multiple curve
design to help
navigate right side
anatomy



Longer braided core
for extended reach

Worley™ Advanced RCSG

Right Sided Coronary Sinus Guide

PROMOTING RELIABLE ACCESS & ACCURATE DELIVERY

To meet the challenges of right sided access during a cardiac resynchronization therapy (CRT), you need an advanced LV delivery system that gives you options to match the tortuosity of your patient's anatomy. The Worley Advanced Right Sided CSG is designed to provide you the options you need to reliably and accurately place your LV leads at the target location.

PRODUCT FEATURES

MULTIPLE SHEATH CURVES to provide back wall support while navigating the right subclavian venous anatomy

BRAIDED CORE with its longer variable shape design provides control and flexibility

DEPTH POSITION MARKER for alignment of the distal tips of the braided core and the RCSG

INTEGRATED PEEL-AWAY HEMOSTASIS VALVE handles like a standard peel-away introducer to reduce risk of blood loss and air embolism


SIDE PORT allows open access infusion and contrast injection during implantation

CRUSH RESISTANT PEBAX®* SHEATH minimizes kinking and facilitates easy removal

9 French ID provides ample room to perform various advanced CRT implantation techniques

SOFT RADIOPAQUE TIP reduces risk of vessel wall trauma and improves fluoroscopic visualisation

VARIABLE SHAPE DESIGN

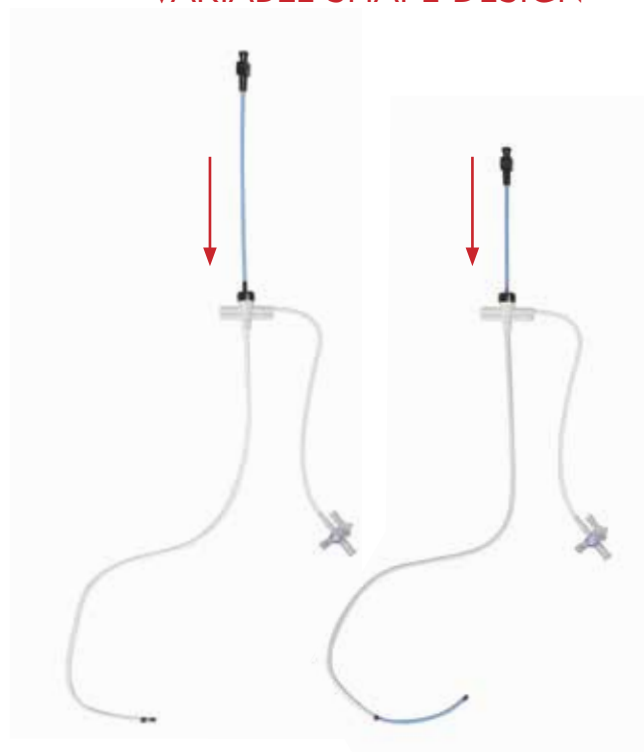

Part Number	French Size	Length
WOR-CSG-BR1-09	9F	40 cm

Kit Contents:

- Splittable Sheath
- Dilator
- Braided Core
- Guidewires:
 - 0.035" x 135 cm
 - 0.032" x 180 cm
- Transvalvular Insertion Tool
- 18G 7 cm Needle
- Syringe

Before using refer to Instructions for Use for indications, contraindications, warnings, precautions, and directions for use.



Understand. Innovate. Deliver.™

merit.com

Merit Medical Systems, Inc.
1600 West Merit Parkway
South Jordan, Utah 84095
1.801.253.1600
1.800.35.MERIT

Merit Medical Europe, Middle
East & Africa (EMEA)
Amerikalaan 42, 6199 AE
Maastricht-Airport
The Netherlands
+31 43 358 82 22

Merit Medical Ireland Ltd.
Parkmore Business Park West
Galway, Ireland
+353 (0) 91 703 733