



Preferential Path. Targeted Delivery.™


Osseoflex® SB

10 gauge compatible Inflatable Bone Tamp

The Osseoflex® SB family of inflatable bone tamps for vertebral augmentation includes both steerable and straight balloons. The Osseoflex SB steerable balloon enables access to the entire vertebral body with unprecedented navigation and precision. The steerable and straight 10 gauge compatible SB balloons are available in 2ml and 4ml balloon sizes.

The Osseoflex SB 10 gauge compatible product line accommodates treatment of compression fractures of vertebra. The Osseoflex SB can be used in either straight or steerable configurations. The SB steerable balloon articulates to follow the channel created by the Osseoflex® SN steerable needle for precision cavity creation.

Utilizing either a unipedicular or bipedicular approach, the Osseoflex SB creates a cavity by compacting cancellous bone for the treatment of vertebral compression fractures.

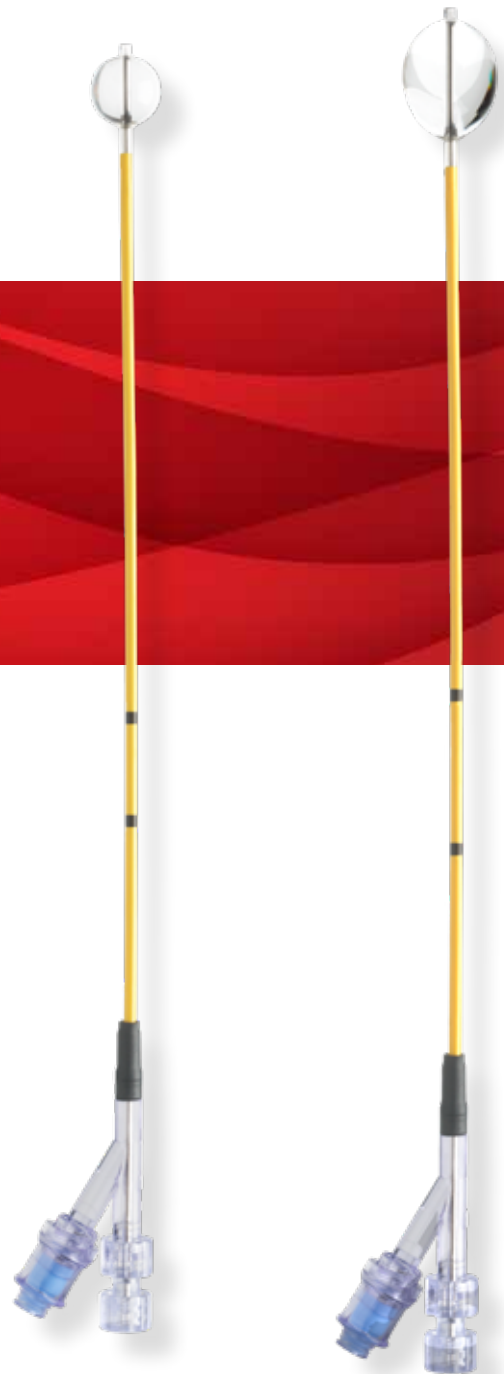

10ga compatible/2ml
Straight Balloon (OT-0222)

10ga compatible/4ml
Straight Balloon (OT-0224)

FEATURES AND BENEFITS

- Radiopaque markers indicate the proximal and distal end of the balloon for accurate positioning prior to inflation*
- Compatible with 10 gauge cannulas
- A distal marker band on the shaft indicates when the balloon tip is flush with the end of the cannula*
- The proximal shaft marker band indicates when the balloon is completely inserted and has cleared the cannula*



*When used with Osseoflex introducer.

IOS

Images may vary, refer to IFU. Images may not be to scale.

Preferential Path. Targeted Delivery.™

Osseoflex SB Steerable Balloons can create a central cavity across the sagittal midline using a unipedicular approach.



OSSEOFLEX SB BALLOONS — 10GA COMPATIBLE

	Part Number	Description	Contents
STEERABLE	OF-0222	Osseoflex SB Steerable Balloon (10ga/2ml)	(1) 10 gauge compatible, 2ml steerable balloon
	OF-0224	Osseoflex SB Steerable Balloon (10ga/4ml)	(1) 10 gauge compatible, 4ml steerable balloon
STRAIGHT	OT-0222	Osseoflex SB Straight Balloon (10ga/2ml)	(1) 10 gauge compatible, 2ml straight balloon
	OT-0224	Osseoflex SB Straight Balloon (10ga/4ml)	(1) 10 gauge compatible, 4ml straight balloon
	IN8VCF	StabiliT Touch™ Syringe	(1) Inflation Device

As with most surgical procedures, there are risks associated with the vertebroplasty, kyphoplasty and vertebral augmentation procedures, including serious complications. For complete information regarding risks, contraindications, warnings, precautions, and adverse events please review the System's Instruction for Use.



Understand. Innovate. Deliver.™

Merit Medical Systems, Inc.
1600 West Merit Parkway
South Jordan, Utah 84095
1.801.208.4300
1.800.35.MERIT

Merit Medical Europe, Middle East, & Africa (EMEA)
Amerikalaan 42, 6199 AE Maastricht-Airport
The Netherlands
+31 43 358 82 22

Merit Medical Ireland Ltd.
Parkmore Business Park West
Galway, Ireland
+353 (0) 91 703 733

merit.com

Australia
+61 (0) 3 8373 4943

Austria
0800 295 374

Belgium
0800 72 906 (Dutch)
0800 73 172 (Français)

Brazil
+55 11 4561 - 0062

Canada
1.800.364.4370

China
+86-10 8561 0788

Denmark
80 88 00 24

Finland
0800 770 586

France
0800 91 60 30

Germany
0800 182 0871

Hong Kong
+852-2207 6888

India
+91-80-41223376

Ireland (Republic)
1800 553 163

Italy
800 897 005

Korea
+8210 2205 0707

Luxembourg
8002 25 22

Netherlands
0800 022 81 84

Norway
800 11629

Portugal
308 801 034

Russia
+7 495 221 89 02

Spain
+34 911238406

Sweden
020 792 445

Switzerland
(Deutsch)
+41 22 518 02 30

(Français)
+41 22 518 02 52

(Italiano)
+41 22 518 00 35

UK
0800 973 115