



**UNIFORMITY.
SIMPLICITY.
VARIETY.**

Fountain® ValveTip

Infusion Systems

Uniform Dispersion • Simplified Procedure • Variety of Options


Fountain® ValveTip

Infusion Systems

UNIFORM AND EXCLUSIVE DISPERSION

- Patented laser-drilled gradient infusion technology
- Facilitates both pulsed-infusion and continuous slow-drip techniques
- 0.035" (0.89mm) guidewire compatible

FOUNTAIN VALVETIP

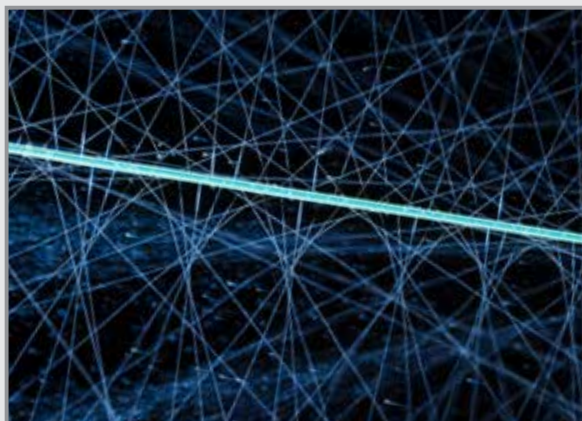

The side holes on the Fountain Infusion Catheter rotate around the catheter to create a more uniform dispersion of therapeutic agent.

COMPETITION



Other competitive catheters have side holes or slits lined up in 4 parallel rows.

SIMPLIFIED PROCEDURE

- ValveTip eliminates the use of an occluding wire
- ValveTip technology minimizes procedural time by eliminating steps associated with an occluding wire
- State-of-the-art valve design streamlines procedure

ACCESSPLUS™ HEMOSTASIS VALVE

Included with every catheter



SQUIRT® FLUID DISPENSING SYSTEM

Compatible



Fountain® ValveTip

Infusion Systems

VARIETY OF OPTIONS

- 4F and 5F configurations
- Wide variety of catheter and infusion lengths
- Variety of options to provide users the products they need

ORDERING INFORMATION

Catalog Number	French Size	Catheter Length	Infusion Segment
IS4-45-10-V	4F	45 cm	10 cm
IS4-45-20-V	4F	45 cm	20 cm
IS4-90-5-V	4F	90 cm	5 cm
IS4-90-10-V	4F	90 cm	10 cm
IS4-90-20-V	4F	90 cm	20 cm
IS4-90-30-V	4F	90 cm	30 cm
IS4-90-40-V	4F	90 cm	40 cm
IS4-90-50-V	4F	90 cm	50 cm
IS4-135-5-V	4F	135 cm	5 cm
IS4-135-10-V	4F	135 cm	10 cm
IS4-135-20-V	4F	135 cm	20 cm
IS4-135-30-V	4F	135 cm	30 cm
IS4-135-40-V	4F	135 cm	40 cm
IS4-135-50-V	4F	135 cm	50 cm
IS5-45-10-V	5F	45 cm	10 cm
IS5-45-20-V	5F	45 cm	20 cm
IS5-90-5-V	5F	90 cm	5 cm
IS5-90-10-V	5F	90 cm	10 cm
IS5-90-20-V	5F	90 cm	20 cm
IS5-90-30-V	5F	90 cm	30 cm
IS5-90-40-V	5F	90 cm	40 cm
IS5-90-50-V	5F	90 cm	50 cm
IS5-135-5-V	5F	135 cm	5 cm
IS5-135-10-V	5F	135 cm	10 cm
IS5-135-20-V	5F	135 cm	20 cm
IS5-135-30-V	5F	135 cm	30 cm
IS5-135-40-V	5F	135 cm	40 cm
IS5-135-50-V	5F	135 cm	50 cm

SQUIRT FLUID DISPENSING SYSTEM

Catalog Number

FDS-100

Before using refer to Instructions for Use for indications, contraindications, warnings, precautions, and directions for use.



Understand. Innovate. Deliver.™

merit.com

Merit Medical Systems, Inc.
1600 West Merit Parkway
South Jordan, Utah 84095
Main: +1.801.253.1600
Customer Service: +1.801.208.4300

Merit Medical Europe, Middle
East, & Africa (EMEA)
Amerikalaan 42, 6199 AE
Maastricht-Airport
The Netherlands
+31 43 358 82 22

Merit Medical Ireland Ltd.
Parkmore Business Park West
Galway, Ireland
+353 (0) 91 703 733

405893001_001 ID 041019