

Cultura™ COLLECTION AND TRANSPORT SYSTEM

The Cultura™ Collection and Transport System is designed with flexibility and convenience to meet your viral testing needs. The system includes a sterile nylon-flocked collection swab with a flexible plastic shaft, a vial with 3.0 mL of viral transport medium, and a resealable biohazard specimen bag with absorbent pad. The sample rapidly and efficiently elutes from the flocked swab when placed in the viral transport medium for testing. The formulation is specially designed for virus preservation and allows for either room temperature or refrigerated transport. Cultura is designed and manufactured in Utah, USA.




THE CULTURA COLLECTION AND TRANSPORT SYSTEM INCLUDES:

CULTURA™ FLEX SWAB (with 100 mm Breakpoint)

CULTURA™ TRANSPORT SYSTEM

BIOHAZARD SPECIMEN BAG (with Absorbent Pad)

Catalog Number	Description	Packaged
VCTS100	1–Sterile flex swab, 100 mm breakpoint 1–Vial with 3 mL viral transport medium 1–Biohazard specimen bag with absorbent pad	100 per box

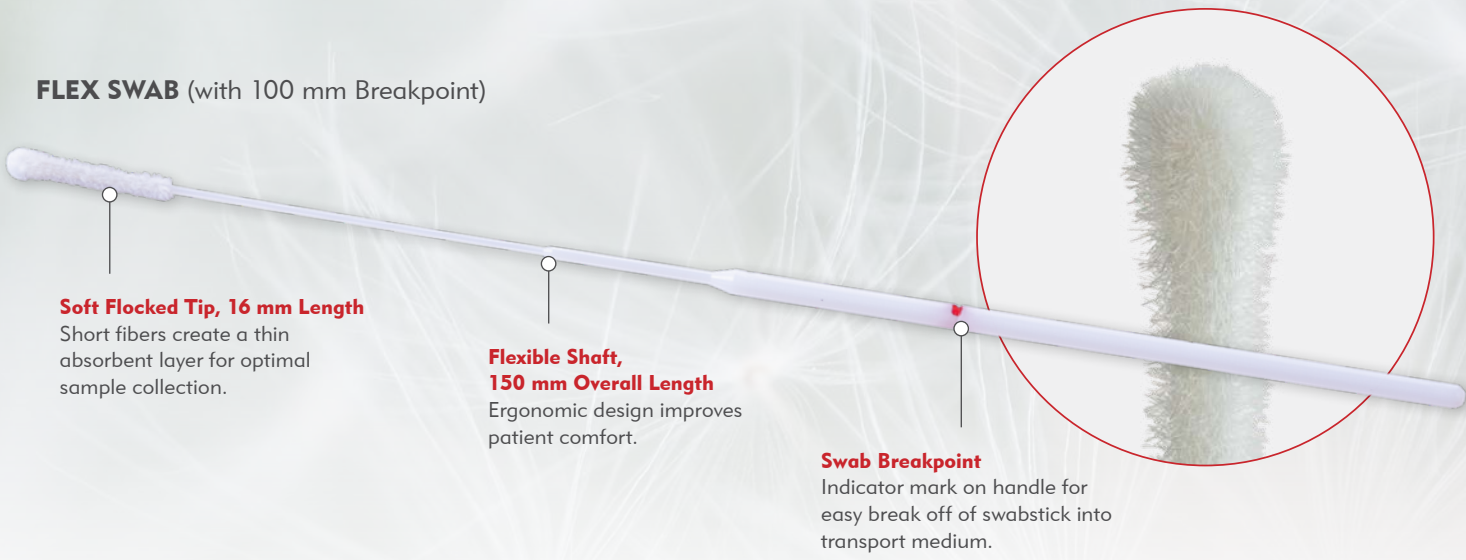

FLEX SWAB (with 100 mm Breakpoint)

Soft Flocked Tip, 16 mm Length

Short fibers create a thin absorbent layer for optimal sample collection.

Flexible Shaft, 150 mm Overall Length

Ergonomic design improves patient comfort.

Swab Breakpoint

Indicator mark on handle for easy break off of swabstick into transport medium.

VIRAL TRANSPORT MEDIUM, 3 mL

- **Temperature Range Flexibility** Room temperature or refrigerated (up to 120 hours) after specimen collection.
- **Maintains Sample Integrity** Formulated with proteins for virus preservation and antibiotics to inhibit bacteria and fungi overgrowth.
- **pH Specific** Phenol red pH indicator ensures medium integrity. Buffer solution maintains neutral pH.
- **Reduces Risk of Biohazard Spills** Self-standing, flat-bottom conical tube.
- **Prevents Leakage** Twist-on cap with dual-sealing surface design.
- **Traceable** Unique QR code on each vial enables efficient tracking of samples in Laboratory Information Management System software.



BIOHAZARD SPECIMEN BAG WITH ABSORBENT PAD, 6"X9"

- **Clean Transport** Absorbent pad allows for specimens to be transported safely and complies with IATA regulations.
- **Simplified Shipping** Outer pouch keeps paperwork separate and dry for shipping.
- **Streamlined Processing** Resealable bag includes pre-printed biohazard symbol with checkboxes for processing.

Before using refer to Instructions for Use for indications, contraindications, warnings, precautions, and directions for use.



Understand. Innovate. Deliver.™

merit.com

Merit Medical Systems, Inc.
1600 West Merit Parkway
South Jordan, Utah 84095
Main: +1.801.253.1600
Customer Service: +1.801.208.4300

Merit Medical Europe, Middle East, & Africa (EMEA)
Amerikalaan 42, 6199 AE
Maastricht-Airport
The Netherlands
+31 43 358 82 22

Merit Medical Ireland Ltd.
Parkmore Business Park West
Galway, Ireland
+353 (0) 91 703 733