



**PUSHABILITY.
TRACKABILITY.
KINK RECOVERY.**

MERIT MAESTRO[®]
Microcatheter

AVAILABLE IN RADIAL ACCESS LENGTH.




PUSHABILITY. TRACKABILITY. KINK RECOVERY.

During your procedures, having a reliable microcatheter will enhance efficiency and provide the outcomes you expect for small-vessel embolization.

RADIAL ACCESS LENGTH

of 165 cm provides the same product benefits as other lengths.

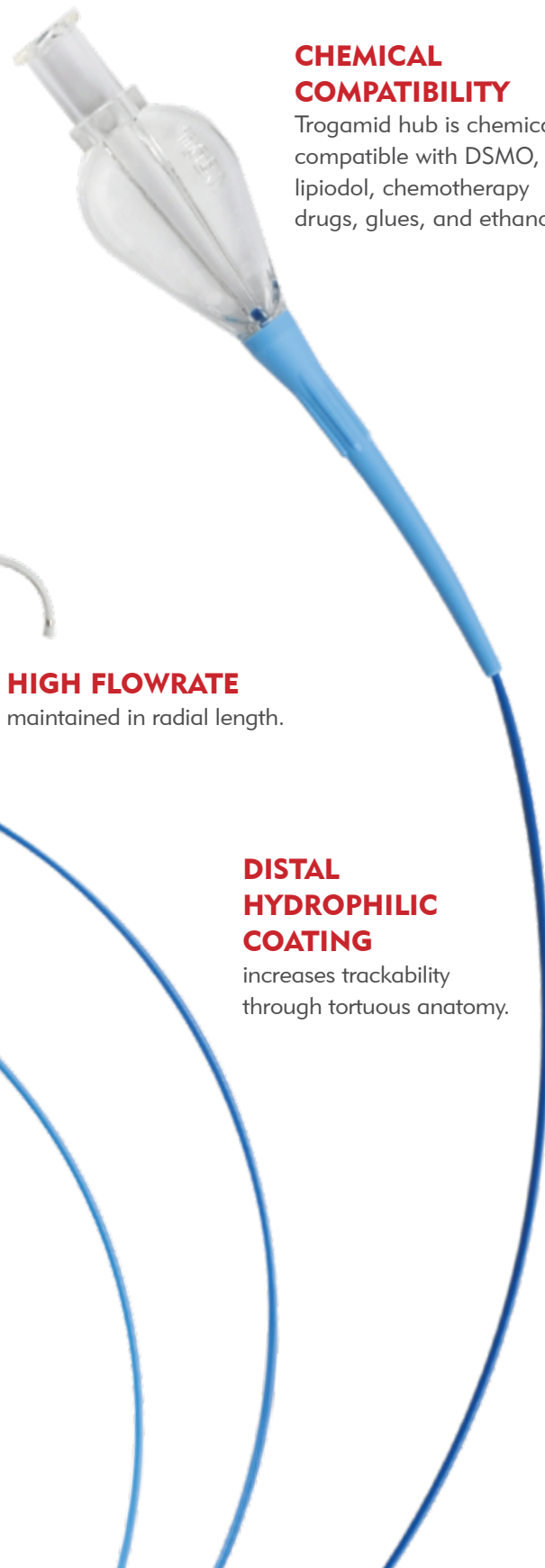

CHEMICAL COMPATIBILITY


Trogamid hub is chemically compatible with DMSO, lipiodol, chemotherapy drugs, glues, and ethanol.

HIGH FLOWRATE

maintained in radial length.

DISTAL HYDROPHILIC COATING

increases trackability through tortuous anatomy.

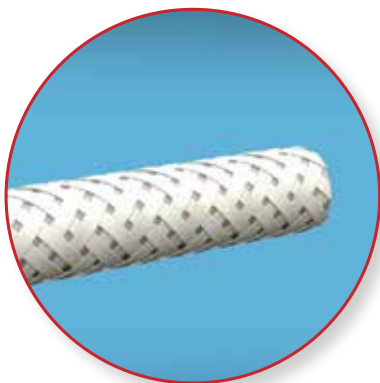

NYLON-RIBBON BRAIDING

provides excellent pushability, trackability, and kink recovery when crossing small tortuous vessels.

Navigating Complex Anatomy

**SOFT, FLEXIBLE
DISTAL TIP**

STRAIGHT

SWAN NECK

45°



COMPLETE SOLUTIONS

Proven for successful delivery of Embosphere® Microspheres.

ORDERING INFORMATION

Lengths with a catalog number in **red** are designed for radial access.

Catalog Number	Maestro Proximal OD	Maestro Distal OD	Maestro Usable Length	Maestro Tip Shape	Maestro Hydrophilic Coating Length	Maestro Distance from Tip to Markerband	Minimum Guide Cath. ID	Maestro ID	Max. Guide Wire Size	Maestro Priming Vol.
----------------	---------------------	-------------------	-----------------------	-------------------	------------------------------------	---	------------------------	------------	----------------------	----------------------

2.8F TAPERED TO 2.1F

28MC21110ST	2.8F	2.1F	110 cm	Straight	80 cm	1.3 mm	0.040" (1.02 mm)	0.018" (0.46 mm)	0.016" (0.41 mm)	0.63 mL
28MC21130ST	2.8F	2.1F	130 cm	Straight	80 cm	1.3 mm	0.040" (1.02 mm)	0.018" (0.46 mm)	0.016" (0.41 mm)	0.70 mL
28MC21150ST	2.8F	2.1F	150 cm	Straight	80 cm	1.3 mm	0.040" (1.02 mm)	0.018" (0.46 mm)	0.016" (0.41 mm)	0.76 mL
28MC21165ST	2.8F	2.1F	165 cm	Straight	80 cm	1.3 mm	0.040" (1.02 mm)	0.018" (0.46 mm)	0.016" (0.41 mm)	0.81 mL
28MC2111045	2.8F	2.1F	110 cm	45 degree	80 cm	1.3 mm	0.040" (1.02 mm)	0.018" (0.46 mm)	0.016" (0.41 mm)	0.63 mL
28MC2113045	2.8F	2.1F	130 cm	45 degree	80 cm	1.3 mm	0.040" (1.02 mm)	0.018" (0.46 mm)	0.016" (0.41 mm)	0.70 mL
28MC2115045	2.8F	2.1F	150 cm	45 degree	80 cm	1.3 mm	0.040" (1.02 mm)	0.018" (0.46 mm)	0.016" (0.41 mm)	0.76 mL
28MC2116545	2.8F	2.1F	165 cm	45 degree	80 cm	1.3 mm	0.040" (1.02 mm)	0.018" (0.46 mm)	0.016" (0.41 mm)	0.81 mL
28MC21110SN	2.8F	2.1F	110 cm	Swan	80 cm	1.3 mm	0.040" (1.02 mm)	0.018" (0.46 mm)	0.016" (0.41 mm)	0.63 mL
28MC21130SN	2.8F	2.1F	130 cm	Swan	80 cm	1.3 mm	0.040" (1.02 mm)	0.018" (0.46 mm)	0.016" (0.41 mm)	0.70 mL
28MC21150SN	2.8F	2.1F	150 cm	Swan	80 cm	1.3 mm	0.040" (1.02 mm)	0.018" (0.46 mm)	0.016" (0.41 mm)	0.76 mL
28MC21165SN	2.8F	2.1F	165 cm	Swan	80 cm	1.3 mm	0.040" (1.02 mm)	0.018" (0.46 mm)	0.016" (0.41 mm)	0.81 mL

2.8F TAPERED TO 2.4F

28MC24110ST	2.8F	2.4F	110 cm	Straight	80 cm	1.3 mm	0.040" (1.02 mm)	0.020" (0.52 mm)	0.018" (0.46 mm)	0.63 mL
28MC24130ST	2.8F	2.4F	130 cm	Straight	80 cm	1.3 mm	0.040" (1.02 mm)	0.020" (0.52 mm)	0.018" (0.46 mm)	0.70 mL
28MC24150ST	2.8F	2.4F	150 cm	Straight	80 cm	1.3 mm	0.040" (1.02 mm)	0.020" (0.52 mm)	0.018" (0.46 mm)	0.76 mL
28MC24165ST	2.8F	2.4F	165 cm	Straight	80 cm	1.3 mm	0.040" (1.02 mm)	0.020" (0.52 mm)	0.018" (0.46 mm)	0.81 mL
28MC2411045	2.8F	2.4F	110 cm	45 degree	80 cm	1.3 mm	0.040" (1.02 mm)	0.020" (0.52 mm)	0.018" (0.46 mm)	0.63 mL
28MC2413045	2.8F	2.4F	130 cm	45 degree	80 cm	1.3 mm	0.040" (1.02 mm)	0.020" (0.52 mm)	0.018" (0.46 mm)	0.70 mL
28MC2415045	2.8F	2.4F	150 cm	45 degree	80 cm	1.3 mm	0.040" (1.02 mm)	0.020" (0.52 mm)	0.018" (0.46 mm)	0.76 mL
28MC2416545	2.8F	2.4F	165 cm	45 degree	80 cm	1.3 mm	0.040" (1.02 mm)	0.020" (0.52 mm)	0.018" (0.46 mm)	0.81 mL
28MC24110SN	2.8F	2.4F	110 cm	Swan	80 cm	1.3 mm	0.040" (1.02 mm)	0.020" (0.52 mm)	0.018" (0.46 mm)	0.63 mL
28MC24130SN	2.8F	2.4F	130 cm	Swan	80 cm	1.3 mm	0.040" (1.02 mm)	0.020" (0.52 mm)	0.018" (0.46 mm)	0.70 mL
28MC24150SN	2.8F	2.4F	150 cm	Swan	80 cm	1.3 mm	0.040" (1.02 mm)	0.020" (0.52 mm)	0.018" (0.46 mm)	0.76 mL
28MC24165SN	2.8F	2.4F	165 cm	Swan	80 cm	1.3 mm	0.040" (1.02 mm)	0.020" (0.52 mm)	0.018" (0.46 mm)	0.81 mL

Catalog Number	Maestro Proximal OD	Maestro Distal OD	Maestro Usable Length	Maestro Tip Shape	Maestro Hydrophilic Coating Length	Maestro Distance from Tip to Markerband	Minimum Guide Cath. ID	Maestro ID	Max. Guide Wire Size	Maestro Priming Vol.
2.8F NO TAPER										
28MC28110ST	2.8F	2.8F	110 cm	Straight	80 cm	1.3 mm	0.040" (1.02 mm)	0.024" (0.62 mm)	0.021" (0.53 mm)	0.63 mL
28MC28130ST	2.8F	2.8F	130 cm	Straight	80 cm	1.3 mm	0.040" (1.02 mm)	0.024" (0.62 mm)	0.021" (0.53 mm)	0.70 mL
28MC28150ST	2.8F	2.8F	150 cm	Straight	80 cm	1.3 mm	0.040" (1.02 mm)	0.024" (0.62 mm)	0.021" (0.53 mm)	0.77 mL
28MC28165ST	2.8F	2.8F	165 cm	Straight	80 cm	1.3 mm	0.040" (1.02 mm)	0.024" (0.62 mm)	0.021" (0.53 mm)	0.82 mL
28MC2811045	2.8F	2.8F	110 cm	45 degree	80 cm	1.3 mm	0.040" (1.02 mm)	0.024" (0.62 mm)	0.021" (0.53 mm)	0.63 mL
28MC2813045	2.8F	2.8F	130 cm	45 degree	80 cm	1.3 mm	0.040" (1.02 mm)	0.024" (0.62 mm)	0.021" (0.53 mm)	0.70 mL
28MC2815045	2.8F	2.8F	150 cm	45 degree	80 cm	1.3 mm	0.040" (1.02 mm)	0.024" (0.62 mm)	0.021" (0.53 mm)	0.77 mL
28MC2816545	2.8F	2.8F	165 cm	45 degree	80 cm	1.3 mm	0.040" (1.02 mm)	0.024" (0.62 mm)	0.021" (0.53 mm)	0.82 mL
28MC28110SN	2.8F	2.8F	110 cm	Swan	80 cm	1.3 mm	0.040" (1.02 mm)	0.024" (0.62 mm)	0.021" (0.53 mm)	0.63 mL
28MC28130SN	2.8F	2.8F	130 cm	Swan	80 cm	1.3 mm	0.040" (1.02 mm)	0.024" (0.62 mm)	0.021" (0.53 mm)	0.70 mL
28MC28150SN	2.8F	2.8F	150 cm	Swan	80 cm	1.3 mm	0.040" (1.02 mm)	0.024" (0.62 mm)	0.021" (0.53 mm)	0.77 mL
28MC28165SN	2.8F	2.8F	165 cm	Swan	80 cm	1.3 mm	0.040" (1.02 mm)	0.024" (0.62 mm)	0.021" (0.53 mm)	0.82 mL

2.9F NO TAPER										
29MC29110ST	2.9F	2.9F	110 cm	Straight	80 cm	1.3 mm	0.042" (1.07 mm)	0.027" (0.68 mm)	0.021" (0.53 mm)	0.69 mL
29MC29130ST	2.9F	2.9F	130 cm	Straight	80 cm	1.3 mm	0.042" (1.07 mm)	0.027" (0.68 mm)	0.021" (0.53 mm)	0.77 mL
29MC29150ST	2.9F	2.9F	150 cm	Straight	80 cm	1.3 mm	0.042" (1.07 mm)	0.027" (0.68 mm)	0.021" (0.53 mm)	0.85 mL
29MC29165ST	2.9F	2.9F	165 cm	Straight	80 cm	1.3 mm	0.042" (1.07 mm)	0.027" (0.68 mm)	0.021" (0.53 mm)	0.90 mL
29MC2911045	2.9F	2.9F	110 cm	45 degree	80 cm	1.3 mm	0.042" (1.07 mm)	0.027" (0.68 mm)	0.021" (0.53 mm)	0.69 mL
29MC2913045	2.9F	2.9F	130 cm	45 degree	80 cm	1.3 mm	0.042" (1.07 mm)	0.027" (0.68 mm)	0.021" (0.53 mm)	0.77 mL
29MC2915045	2.9F	2.9F	150 cm	45 degree	80 cm	1.3 mm	0.042" (1.07 mm)	0.027" (0.68 mm)	0.021" (0.53 mm)	0.85 mL
29MC2916545	2.9F	2.9F	165 cm	45 degree	80 cm	1.3 mm	0.042" (1.07 mm)	0.027" (0.68 mm)	0.021" (0.53 mm)	0.90 mL
29MC29110SN	2.9F	2.9F	110 cm	Swan	80 cm	1.3 mm	0.042" (1.07 mm)	0.027" (0.68 mm)	0.021" (0.53 mm)	0.69 mL
29MC29130SN	2.9F	2.9F	130 cm	Swan	80 cm	1.3 mm	0.042" (1.07 mm)	0.027" (0.68 mm)	0.021" (0.53 mm)	0.77 mL
29MC29150SN	2.9F	2.9F	150 cm	Swan	80 cm	1.3 mm	0.042" (1.07 mm)	0.027" (0.68 mm)	0.021" (0.53 mm)	0.85 mL
29MC29165SN	2.9F	2.9F	165 cm	Swan	80 cm	1.3 mm	0.042" (1.07 mm)	0.027" (0.68 mm)	0.021" (0.53 mm)	0.90 mL

Packaged 1 unit per box.

CATHETER COMPATIBILITY

Size	Maximum Microsphere Size	Maximum Particle Size	Maximum Coil Size	Maximum Injection Pressure
2.8F–2.1F	≤ 700 μm	≤ 500 μm	≤ 0.016" (0.41 mm)	800 psi (5,515 kPa)
2.8F–2.4F	≤ 700 μm	≤ 700 μm	≤ 0.018" (0.46 mm)	800 psi (5,515 kPa)
2.8F–2.8F	≤ 700 μm	≤ 700 μm	≤ 0.018" (0.46 mm)	800 psi (5,515 kPa)
2.9F–2.9F	≤ 900 μm	≤ 900 μm	— — — — —	800 psi (5,515 kPa)

Data on File

Before using refer to Instructions for Use for indications, contraindications, warnings, precautions, and directions for use.



Understand. Innovate. Deliver.™

merit.com

Merit Medical Systems, Inc.
1600 West Merit Parkway
South Jordan, Utah 84095
1.801.253.1600
1.800.35.MERIT

Merit Medical Europe, Middle
East & Africa (EMEA)
Amerikalaan 42, 6199 AE
Maastricht-Airport
The Netherlands
+31 43 358 82 22

Merit Medical Ireland Ltd.
Parkmore Business Park West
Galway, Ireland
+353 (0) 91 703 733